



Technology for Vacuum Systems

DIAPHRAGM PUMP

MV 10 VARIO select

MD 12 VARIO select

ME 16 VARIO select



Instructions for use



**Original instructions
Keep for further use!**

This manual is only to be used and distributed in its complete and original form. It is strictly the user's responsibility to carefully check the validity of this manual with respect to the product.

Manufacturer:
VACUUBRAND GMBH + CO KG
Alfred-Zippe-Str. 4
97877 Wertheim
GERMANY

Phone:

- Head office +49 9342 808-0
- Sales +49 9342 808-5550
- Service +49 9342 808-5660

Fax: +49 9342 808-5555
Email: info@vacuubrand.com
Web: www.vacuubrand.com

*Thank you for purchasing this product from **VACUUBRAND GMBH + CO KG** . You have chosen a modern and technically high quality product.*

TABLE OF CONTENT

1	Introduction	5
1.1	User information	5
1.2	About this document.	6
1.2.1	Manual structure	6
1.2.2	Display conventions	7
1.2.3	Symbols and icons	8
1.2.4	Handling instructions (action steps).	9
1.2.5	Abbreviations	9
1.2.6	Term definitions.	10
2	Safety information	11
2.1	Usage.	11
2.1.1	Intended use	11
2.1.2	Improper use.	12
2.1.3	Foreseeable misuse.	12
2.2	Obligations	13
2.2.1	Operator obligations	13
2.2.2	Personnel obligations	13
2.3	Target group description.	14
2.4	General safety information	15
2.4.1	Protective clothing.	15
2.4.2	Safety precautions	15
2.4.3	Laboratory and working materials.	16
2.4.4	Eliminate sources of danger	17
2.5	Motor protection	20
2.6	ATEX equipment category	21
2.7	Proper disposal	22
3	Product description	23
3.1	Schematic design.	23
3.2	Diaphragm pump series	24
4	Installation and connection	25
4.1	Transport	25
4.2	Installation	26
4.3	Controller base	28
4.4	Connection	29
4.4.1	Assemble silencer (EX)	29
4.4.2	Vacuum connection (IN)	30
4.4.3	Outlet connection (EX)	32
4.4.4	Connect venting valve (optional).	34
4.4.5	Electrical connection.	35

5	Commissioning (operation)	37
5.1	Switch on	37
5.2	Operation	37
5.2.1	Operation (→ see description of controller).	39
5.3	Shutdown (switch off).	40
5.4	Storage.	41
6	Troubleshooting	43
6.1	Technical support.	43
6.2	Error – Cause – Remedy	43
7	Cleaning and maintenance	47
7.1	Information on service work	48
7.2	Cleaning	50
7.2.1	Diaphragm pump	50
7.2.2	Clean or replace molded PTFE hoses	50
7.2.3	Clean or replace the controller	50
7.3	Diaphragm pump maintenance	51
7.3.1	Maintenance items	51
7.3.2	Preparation	52
7.3.3	Replacing the diaphragms and valves	55
8	Appendix	70
8.1	Technical information	70
8.1.1	Technical data.	70
8.1.2	Wetted materials.	73
8.1.3	Rating plate.	74
8.2	Ordering information.	75
8.3	Service	77
8.4	Index.	78
8.5	EC Declaration of Conformity	80

1 Introduction

This manual is part of your product. The manual applies to all variants of the pump and is intended in particular for laboratory staff.

1.1 User information

Safety

Instructions for use
and safety

- Read this manual thoroughly and completely before using the product.
- Keep this manual in an easily accessible location.
- Correct use of the product is essential for safe operation. Comply with all safety information provided!
- In addition to this manual, adhere to the accident prevention regulations and industrial safety regulations applicable in the country of use.

General

General
information

- For easier readability, the general term *diaphragm pump is used as an equivalent to and instead of the product name Mx 1x VARIO select diaphragm pump.*
- If passing the product on to a third party, also give them this manual.
- The illustrations in this manual are only intended to facilitate comprehension.
- We reserve the right to make technical changes in the course of continuous product improvement.

Copyright

Copyright ©

The content of this manual is protected by copyright. Only copies for internal use are allowed, e.g., for professional training.

© **VACUUBRAND GMBH + CO KG**

Contact

Contact
US

- If your manual is incomplete, you can request a replacement. Alternatively, you can use our download portal: www.vacuubrand.com
- You are welcome to contact us at any time in writing or by telephone if you would like more information, have questions about our products or wish to share feedback with us.
- When contacting our Service Department, please have the serial number and product type at hand → see *Rating plate on the product*.

1.2 About this document


1.2.1 Manual structure

Modular instructions
for use

The manual has a modular structure with separate instruction modules for the diaphragm pump, vacuum controller, and any accessories.

Instruction modules

Pump series and
instructions for use

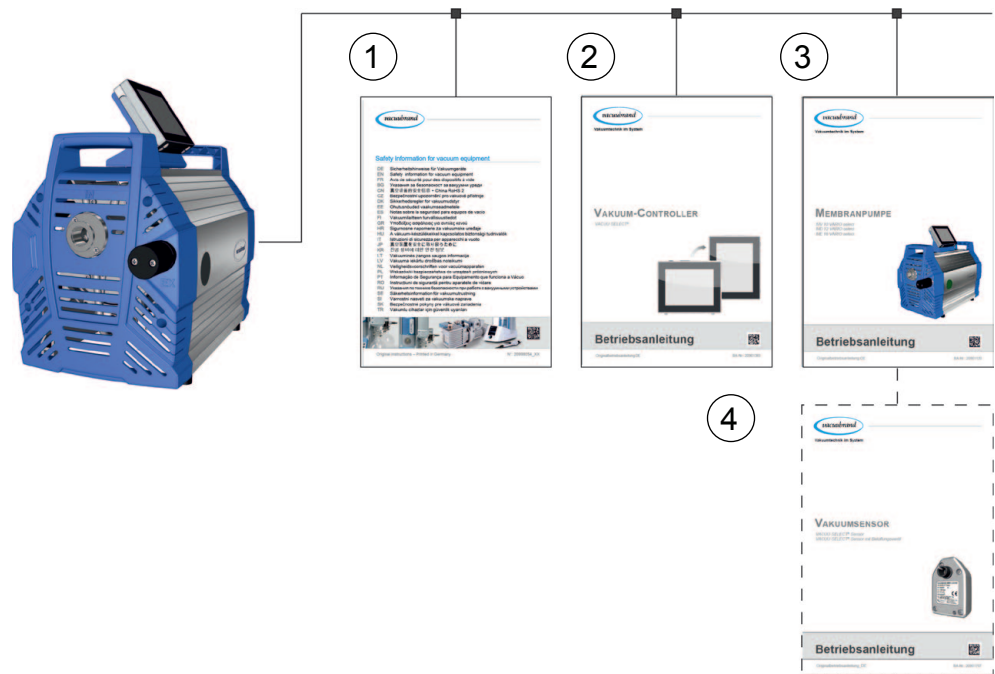



- | | |
|---|--|
| 1 | Safety information for vacuum equipment |
| 2 | Description: Vacuum controller – control and operation |
| 3 | Description: Vacuum pump – connection, operation, maintenance, mechanics |
| 4 | Optional description: Accessories |


1.2.2 Display conventions

Warning levels

Display conventions

	DANGER
	<p>Indicates an imminent hazardous situation. Disregarding the situation will result in serious and even fatal injury or death.</p> <p>⇒ Take appropriate action to avoid dangerous situations!</p>

	WARNING
	<p>Indicates a potentially hazardous situation. Disregarding the situation could result in serious, even fatal injury or massive damage to property.</p> <p>⇒ Take appropriate action to avoid dangerous situations!</p>


	CAUTION
	<p>Indicates a potentially hazardous situation. Disregarding the situation could result in slight or minor injury or damage to property.</p> <p>⇒ Take appropriate action to avoid dangerous situations!</p>

NOTICE	
<p>Indicates a potentially harmful situation. Disregarding the notice could lead to material damage.</p>	

Additional notes

IMPORTANT!

- ⇒ Information or specific recommendation which must be observed.
- ⇒ Important information for proper operation.

	<ul style="list-style-type: none"> ⇒ Helpful tips + tricks ⇒ Additional notes
---	---

1.2.3 Symbols and icons

This manual uses symbols and icons. Safety symbols indicate specific risks associated with handling the product. Symbols and icons are designed to help you identify risks more easily.









Safety symbols

Explanation of safety symbols

	Hazardous substance - hazards to human health.		General prohibition sign.
	General warning sign.		Warning Risk of explosion.
	Danger: electricity.		Danger: hot surface
	General mandatory sign.		Disconnect power plug.
	Wear chemical resistant protective gloves.		Wear protective goggles.
	Read the instructions for repair.		

Additional icons

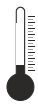Additional symbols

	Positive example – Do this! Result – OK		Negative example – Do not do this!
	Refers to content in this manual.		Refers to content of other supplementary documents.
	Installation at temperatures < 40 °C.		Ensure sufficient air circulation.
	Flow arrow Inlet – vacuum connection		Flow arrow Outlet – exhaust gas

1.2.4 Handling instructions (action steps)

Display of operating steps

Instructions (single step)

⇒ Perform the step described.

- Result of action

Instructions (multiple steps)

1. First step


2. Next step

- Result of action

Perform the steps in the order described.

1.2.5 Abbreviations

Abbreviations

abs.	absolute
ATM	Atmospheric pressure (pressure graphic, program)
d_i (di)	Interior diameter
DN	Nominal diameter
ECTFE	Ethylene chlorotrifluoroethylene
e.g.	for example
ETFE	Ethylene/Tetrafluoroethylene
EX*	Outlet (exhaust, exit), exhaust gas connection
	ATEX equipment labeling
Fig.	Figure
FPM	Fluoroelastomer
IN*	Inlet, vacuum connection
KF	Small flange
max.	Maximum value
min.	Minimum value
PP	Polypropylene
PPS	Polyphenylene sulphide
PTFE	Polytetrafluorethylene
resp.	responsible (supervising)
RMA-N°	Return Merchandise Authorization number
SW	Wrench size (tool)

* Labeling on vacuum pump or component

1.2.6 Term definitions

Product-specific terms

Fine vacuum.	Pressure range in vacuum technology, from: 1 mbar–0,001 mbar
Rough vacuum	Pressure range in vacuum technology, from: Atmospheric pressure–1 mbar
Mx 1x VARIO select	Vacuum pump with variable speed motor for precise vacuum control including VACUU·SELECT® controller and VACUU·SELECT® Sensor .
VACUU·BUS®	Bus system from VACUUBRAND for communication between peripheral devices with VACUU·BUS® enabled gauges and controllers. The maximum permissible cable length is 30 m.
VACUU·BUS® address	Address which enables the VACUU·BUS® client to be unambiguously assigned within the bus system, e.g., for connecting multiple sensors with the same measurement range.
VACUU·BUS® client	Peripheral device or component with VACUU·BUS® port which is integrated in the bus system, e.g., sensors, valves, level indicators, etc.
VACUU·BUS® connector	4-pin round connector for the bus system from VACUUBRAND .
VACUU·BUS® configuration	Assigning a different VACUU·BUS® address to a VACUU·BUS® component using a gauge or controller.
VACUU·LAN®	Local area vacuum network.
VACUU·SELECT®	Vacuum controller, controller with touchscreen; consisting of operating panel and vacuum sensor.
VACUU·SELECT® Sensor	External vacuum sensor ▶ for VACUU·SELECT® <i>or</i> ▶ separately as an independent vacuum sensor.
VARIO® drive	Speed control for vacuum pump; the motor runs only as fast as necessary to meet demand.

2 Safety information

The information in this chapter must be observed by everyone who works with the product described here.

The safety instructions are valid for the complete life cycle of the product.

2.1 Usage

Only use the product if it is in perfect working condition.

2.1.1 Intended use

Intended use A diaphragm pump of the *Mx 1x VARIO select* product series is a vacuum system consisting of a vacuum pump, controller and sensor to create and control rough vacuum in designated systems, e.g., as backing pump for high vacuum pumps, for vacuum drying or in systems with VACUU-LAN local area vacuum network etc. The vacuum system may only be used indoors in a non-explosive atmosphere.

Intended use also includes:



- observing the information in the document **Safety information for vacuum equipment**,
- observing the manual,
- observing the manual of connected components,
- observing the inspection and maintenance intervals and having maintenance performed by appropriately qualified personnel.
- using only approved accessories or spare parts.

Any other use is considered improper use.

2.1.2 Improper use

Improper use Incorrect use or any application which does not correspond to the technical data may result in injury or damage to property.

Improper use includes:

- using the product contrary to its intended use,
- using the product in non-commercial environments, unless the necessary protective measures and precautions have been taken by the company,
- operation under inadmissible environmental and operating conditions,
- operation despite obvious faults or defective safety devices,
- unauthorized extensions or conversions, in particular when these impair safety,
- usage despite incomplete assembly,
- operation with sharp-edged objects,
- pulling plug-in connections on the cable out of the socket,
- aspirating, conveying, or compressing solids or fluids.

2.1.3 Foreseeable misuse

Foreseeable misuse In addition to improper use, there are types of use which are prohibited when handling the product:

Prohibited types of use are, in particular:



- use on humans or animals,
- installation and operation in potentially explosive atmospheres,
- use in mines or underground,
- using the product to generate pressure,
- fully exposing vacuum equipment to the vacuum,
- immersing vacuum equipment in liquids, or exposing it to water spray or steam jets,

Foreseeable misuse

- pumping oxidizing and pyrophoric substances, liquids or solids,
- pumping hot, unstable, or explosive media,
- pumping substances which may react explosively under impact and/or elevated temperature without an air supply.

IMPORTANT!

No foreign bodies, hot gases or flames from the application must be allowed to enter the equipment.

2.2 Obligations

2.2.1 Operator obligations

Operator obligations

The owner defines the responsibilities and ensures that only trained personnel or specialists work at the vacuum system. This applies in particular to connection, assembly and maintenance work, and troubleshooting.

Users in the areas of competence in the *Responsibility matrix* must possess the relevant qualifications for the activities listed. In particular work on electrical equipment must be performed only by qualified electricians.

2.2.2 Personnel obligations

Personnel obligations

In the case of activities which require protective clothing, personal protective equipment as specified by the operator is to be worn. If the vacuum system is not in proper working order, it must be prevented from being accidentally switched back on.

- ⇒ Always be conscious of safety and work in a safe manner.
- ⇒ Observe instructions issued by the operator, and national regulations on accident prevention and industrial safety.



The way individuals act can help to prevent accidents at work.

2.3 Target group description

Target groups The manual must be read and observed by every person who is tasked with the activities described below.

Personnel qualification

Qualification description

Operator	Laboratory staff, such as chemists, laboratory technicians
Specialist	Person with professional qualification in mechanics, electrical equipment or laboratory devices
Responsible specialist	Similar to a specialist, with additional specialist responsibility, or responsibility for a department or division

Responsibility matrix

Responsibility Assignment Matrix

Task (Job)	Operator	Specialist	Responsible specialist
Installation	x	x	x
Initial use	x	x	x
Network integration			x
Operation	x	x	x
Error report	x	x	x
Remedy	(x)	x	x
Maintenance		x	x
Repair ¹		x	x
Repair order			x
Cleaning, simple	x	x	x
Shutdown	x	x	x
Decontamination ²		x	x

1 see also website:

VACUUBRAND > Support > [Instructions for repair](#)

2 Alternatively, arrange for decontamination by a qualified service provider


2.4 General safety information

Quality standards
and
safety

Products from **VACUUBRAND GMBH + CO KG** are subject to stringent quality testing with regard to safety and operation. Each product undergoes a comprehensive test program prior to delivery.

2.4.1 Protective clothing

Protective clothing

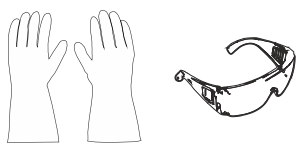

No special protective clothing is required to operate the vacuum pump. Observe instructions issued by the operator for your workplace.

During cleaning, maintenance and repair work, we recommend wearing chemical-resistant protective gloves, protective clothing and protective goggles.

IMPORTANT!

⇒ When handling chemicals, wear your personal protective equipment.

2.4.2 Safety precautions

Safety precautions


- ⇒ Use the vacuum device only if you have understood its function and this manual.
- ⇒ Replace defective parts immediately, e.g., a broken cable, faulty hoses, etc.
- ⇒ Use only original accessories and components which are designed for the vacuum technology, such as a vacuum hose, separator, vacuum valve, etc.
- ⇒ When handling contaminated parts, follow the relevant regulations and safety precautions, this also applies to equipment sent in for repair.

IMPORTANT!

Prior to returning any product to our Service Department for repair, contamination from hazardous substances needs to be excluded.

⇒ Fill out the [Health and Safety Clearance form](#) in full and confirm with your signature.

2.4.3 Laboratory and working materials

	DANGER
	<p>Hazardous substances could be discharged at the outlet.</p> <p>During aspiration, hazardous, toxic substances at the outlet can get into the ambient air.</p> <ul style="list-style-type: none"> ⇒ Observe the national regulations for safe handling of hazardous substances. ⇒ Please note that residual process media may pose a danger to people and the environment. ⇒ Mount and use suitable separators, filters or fume hood devices.

Hazards due to different substances

Pumping different substances

Pumping different substances or media can cause the substances to react with one another.

Working materials which get into the vacuum pump with the gas flow can damage the vacuum pump. Hazardous substances can deposit in the vacuum pump.

Possible protective measures, depending on the application:

- ⇒ Flush the vacuum pump with inert gas or air before changing the medium to be pumped.
- ⇒ Use inert gas to dilute critical mixtures.
- ⇒ Prevent the release of hazardous, toxic, explosive, corrosive fluids, gases or vapors or those that are harmful to health or the environment, for example, through suitable laboratory facilities with a fume hood and ventilation control.
- ⇒ Protect the inside of the vacuum pump from deposits or moisture.
- ⇒ Be aware of interactions and possible chemical reactions of the pumped media.
- ⇒ Check the compatibility of the pumped substances with the wetted materials of the vacuum pump.
- ⇒ Contact us if you have concerns about using your vacuum pump with certain working materials or media.

2.4.4 Eliminate sources of danger

Take mechanical stability into account

Note
mechanical load
capacity

The high compression ratio of the pump may result in a higher pressure at the outlet than the mechanical stability of the system allows.

- ⇒ Always ensure that the outlet line is clear and non-pressurized. The outlet must not be blocked to ensure that gases can exit freely.
- ⇒ Prevent uncontrolled overpressure, e.g., due to a locked or blocked piping system, condensate or clogged outlet line or silencer.
- ⇒ High gas flow can lead to overpressure at the silencer. In case of permanently high gas flow replace the silencer at the outlet by a small flange connection or a hose nozzle and connect an outlet line.
- ⇒ At the gas connections, the connections for the inlet **IN** and outlet **EX** must not be mixed up.
- ⇒ Be aware of the max. pressures at the inlet and outlet of the pump as well as the max. admissible differential pressure between the inlet and outlet, according to. **8.1.1 Technical data on page 70**
- ⇒ The system to be evacuated as well as all hose connections must be mechanically stable.

Prevent condensate return

Prevent backup in
the outlet line

Condensate can damage the pump head. Condensate must not flow back into the outlet **EX** or pump head through the hose line. No liquid should accumulate inside the exhaust hose or inside the silencer.

- ⇒ Avoid condensate return by using a separator (accessory). Condensate must not enter the inside of the housing via the hose lines.
- ⇒ Preferably route the exhaust gas hose with a fall from the outlet, i.e., running downward so that no backup forms.

Prevent foreign bodies inside the pump

Observe vacuum pump dimensioning

Particles, liquids and dust must not get inside the vacuum pump.

- ⇒ Do not pump any substances which could form deposits inside the vacuum pump.
- ⇒ Install suitable separators and/or filters upstream of the inlet. Suitable filters are chemically resistant, clog-proof and have a reliable flow rate, for example.
- ⇒ Replace porous vacuum hoses without delay.

Hazards during venting

Hazards when venting

Depending on the application, explosive mixtures can form or other hazardous situations can arise in systems.

Hazards due to residual energy

Possible residual energies

After the vacuum pump has been switched off and disconnected from the power supply, there may still be dangers due to residual energy:

- Thermal energy: Motor waste heat, hot surface, compression heat.
 - ⇒ Allow the vacuum pump to cool down.
- Electrical energy: The capacitors on the electronic assembly have a discharge time of up to 3 minutes.
 - ⇒ Wait until the capacitors have discharged.

Risk of burns due to hot surfaces or overheating

Surface temperatures

The surface of the vacuum pump can reach operating temperatures > **70 °C**, in particular when pumping heated media. The surface temperature of the silencer in particular might be elevated in case of high gas flow.

- ⇒ Avoid direct contact with the surface.
- ⇒ Use protection against accidental contact if the surface temperature is regularly elevated.
- ⇒ Allow the vacuum pump to cool down before performing maintenance work.

Overheating

The vacuum pump can be damaged due to overheating. Possible causes include insufficient air supply to the fan and failure to maintain minimum distances.

- ⇒ When installing the device, ensure that there is a minimum distance of 5 cm between the fan and adjacent parts (such as the housing, walls, etc.).
- ⇒ Always ensure a sufficient air supply; if applicable, provide external forced ventilation.
- ⇒ Place the device on a stable surface; a soft surface such as foam rubber as a sound absorber can impair and block the air supply.
- ⇒ Clean polluted ventilation slots.
- ⇒ Remove covers from the device before operating it.
- ⇒ Avoid excessive heat input due to hot process gases.
- ⇒ Observe the maximum admissible media temperature
→ *see chapter: 8.1.1 Technical data on page 70.*

Keep signs legible

Signs and labels Keep labels and information symbols and warning labels always in a well readable condition:

- ⇒ Connection labels
- ⇒ Warning signs and notice labels
- ⇒ Motor data and rating plates

2.5 Motor protection

Overheating protection, blockage protection The pump motor has a temperature sensor on the circuit board as overload protection. In the event of excessive temperature or if the motor is blocked, the vacuum pump switches off.

Procedure for switching vacuum pump back on If the vacuum pump is switched off due to these safety precautions, the error must be cleared manually: Unplug pump from the power supply → Eliminate cause of error → Switch vacuum pump back on

2.6 ATEX equipment category

Installation and potentially explosive atmospheres



Installation and operation in areas where potentially explosive atmospheres can develop to a hazardous degree is not permitted.

ATEX approval only applies to the internal, wetted parts of the of the product, not to its surroundings.

ATEX equipment labeling

ATEX equipment category

Vacuum equipment labeled with Ex has ATEX approval in line with the ATEX marking on the rating plate.



- ⇒ Only use the product if it is in perfect working condition.
- ⇒ The devices are designed for a low level of mechanical stress and must be installed in such a way that they cannot sustain mechanical damage from the outside.

ATEX equipment category and peripherals

The ATEX equipment category of the product is dependent on the connected components and peripheral devices. Components and connected peripherals need to have the same or higher ATEX approval.

Prevent ignition sources

The use of venting valves is only permitted if this would not normally, or only rarely, cause explosive mixtures within the device, or do so only for a short time.

- ⇒ If necessary vent with inert gas.


Information on the ATEX equipment category is also available on our website at: www.vacuubrand.com/.../Information-ATEX

2.7 Proper disposal

NOTICE

Risk of environmental damage due to incorrect disposal of the product.

- ⇒ Do not dispose your product in household waste!
Electronic components are subject to hazardous waste treatment and must only be disposed of by certified specialists.
- ⇒ Observe the national regulations for safe disposal and environmental protection.
- ⇒ Detailed information on the respective regulations can be obtained from your local administrative authority.

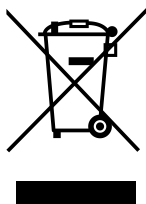


3 Product description

The diaphragm pumps described essentially consist of a diaphragm pump with VARIO® drive, a **VACUU-SELECT®** vacuum controller and a **VACUU-SELECT® Sensor**

3.1 Schematic design

Side and front view

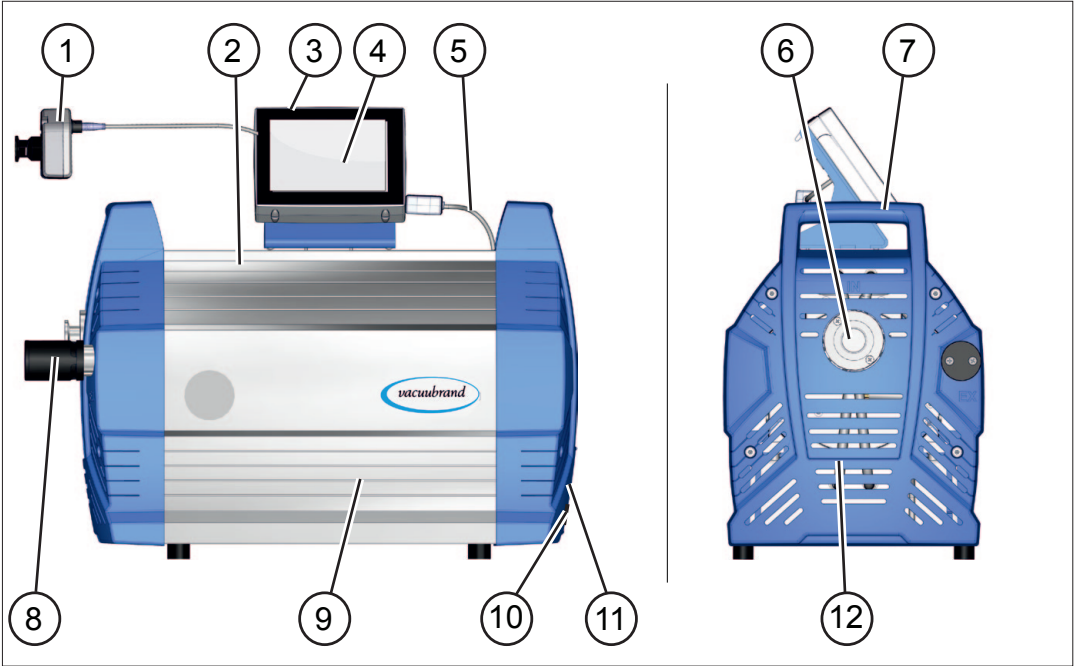

Meaning

1	VACUU-SELECT® Sensor , to be mounted externally on suction line
2	Diaphragm pump
3	Vacuum controller ON/OFF button
4	VACUU-SELECT® operating panel
5	VACUU-SELECT® VACUU-BUS cable (power supply + control cable)
6	Vacuum connection – inlet IN
7	Handles (2x)
8	Silencer– outlet connection – outlet EX
9	Side panel, cover
10	Power connection, ON/OFF button (rocker switch)
11	Rating plate
12	Housing section with handle, front

3.2 Diaphragm pump series

The diaphragm pumps do not differ in their outward appearance.
 → see figure: *3.1 Schematic design on page 23*

The diaphragm pumps differ in the internal connection of the pump heads.

Diaphragm pumps Mx 1x VARIO select

Stages of
diaphragm pump

Diaphragm pump	Pump heads	Stages
ME 16 VARIO select	8	1
MD 12 VARIO select	8	3
MV 10 VARIO select	8	4

4 Installation and connection

4.1 Transport

Products from **VACUUBRAND** are packed in sturdy, recyclable packaging.



The original packaging is accurately matched to your product for safe transport.

⇒ If possible, please keep the original packaging, e. g., for returning the product for repair.


Goods arrival

Check the shipment for transport damage and completeness.

⇒ Immediately report any transport damage in writing to the supplier.

Unpacking

→ Example
Diaphragm pump in
original packaging
with enclosed
packages

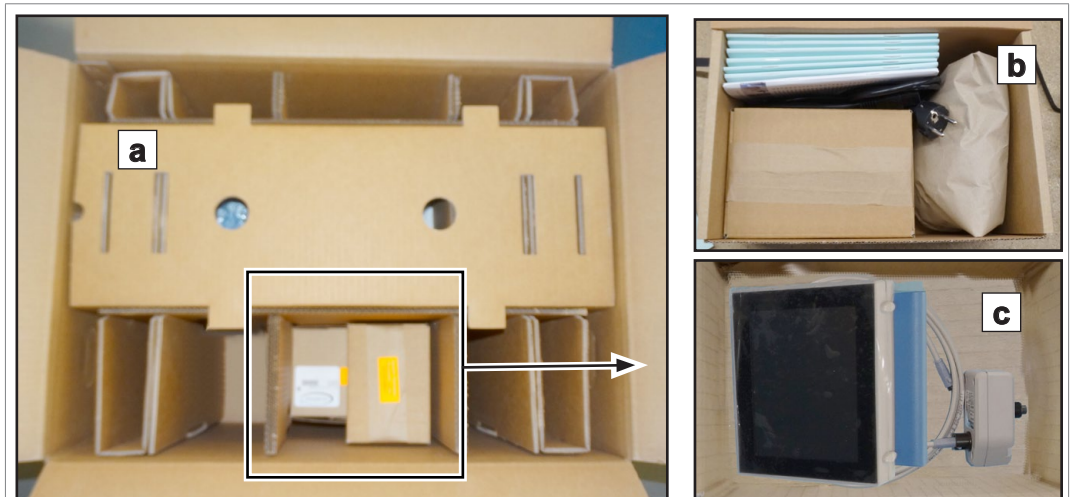

(a) = diaphragm pump
(b) = manual , cable, silencer, and any accessories
(c) = controller, vacuum sensor, cable

⇒ Remove all enclosed packages from their original packaging and unpack them.

⇒ Compare the scope of delivery with the delivery note.

→ Example
Lift out the
diaphragm pump



- ⇒ Please note that a **diaphragm pump can weigh approx. 28 kg**. We recommend using a lifting aid.
- ⇒ Lift the device out of the packaging by the side handles.

4.2 Installation

NOTICE

Condensate can damage the electronics.

A large temperature difference between the storage location and the installation location can cause condensation.

- ⇒ After goods receipt or storage, allow your vacuum device to acclimatize for at least 3-4 hours before initial use.

Check installation conditions


Check installation
conditions

- The device is acclimatized.
- Ambient conditions have been observed and are within the limitation of use.
- The pump must have a stable and secure base without additional mechanical contact apart from the pump feet.

Installing the vacuum pump

⇒ Place the vacuum pump on a stable, non-vibrating, level surface.

→ Example Sketch of minimum distances in laboratory furniture

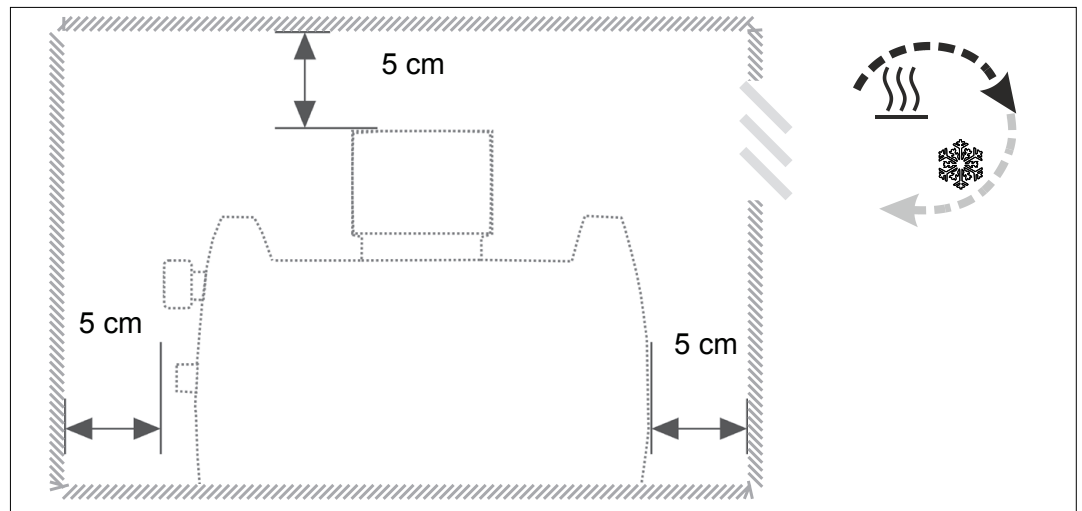

IMPORTANT!

⇒ When installing in lab furniture, maintain a minimum distance of 5 cm (2 in.) to adjacent objects or surfaces.

⇒ Prevent heat accumulation and ensure sufficient air circulation, especially in closed housings.

Observe limitations of use

Observe limitation of use

Limitation of use		(US)
Ambient temperature	10–40 °C	50–104°F
Max. altitude	2000 m above NHN	6562 ft above sea level
Minimum distance to adjacent parts	5 cm	2 in
Relative humidity	30–85 %, non condensing	
Protection class	IP 40/IK 08	
Prevent condensation or contamination from dust, liquids, or corrosive gases.		

IMPORTANT!

⇒ Note the IP protection class. IP protection is only guaranteed if the device is appropriately mounted and connected.

⇒ For connection also note the rating plate data and chapter **8.1.1 Technical data on page 70.**


4.3 Controller base

The base, controller, screw fittings and vacuum sensor are enclosed separately. Before installation, the base can be mounted on the pump and the controller clipped into place.

Alternatively, the controller can be clipped into a recess in the lab furniture or used as a freestanding unit (unfold the stand).

Mount the base


Mount the base to the diaphragm pump (option)

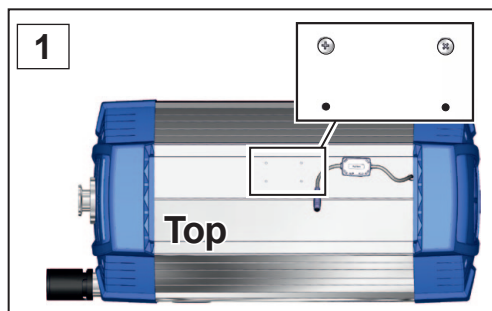

1. Unscrew the screws; Phillips screwdriver size 1.



2. Position the base on the diaphragm pump.

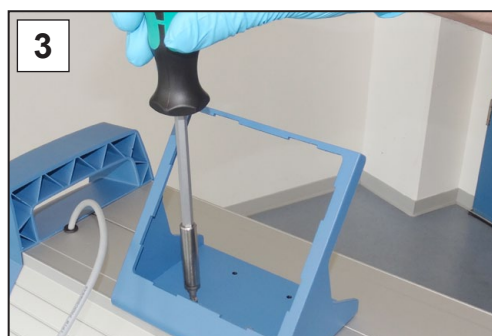

3. Screw the base onto the diaphragm pump.



4. Clip the controller into the base.



5. Plug the VACUU·BUS cable into the power connection on the back of the controller.



6. Also plug in the VACUU·BUS cables from peripheral devices. Use Y adapters (accessories) if there are not enough connections.

4.4 Connection

The diaphragm pumps have a vacuum connection and an outlet connection. Connect your diaphragm pump as described in the examples below.

4.4.1 Assemble silencer (EX)

Silencer at the outlet EX

As standard the diaphragm pumps are equipped with a silencer at the outlet (EX). The silencer is separately packed.



WARNING

Risk of bursting due to internal overpressure at the silencer.

Inadmissibly high pressure at the silencer can cause the vacuum pump to burst or damage the pump bearings, diaphragms and valves.

Internal overpressure may build up in case of high gas flow rate or in case of deposits inside the silencer caused by pumping gases containing dust or solvent vapors.

- ⇒ Do not pump any substances which could form deposits inside the silencer.
- ⇒ In case of permanently high gas flow or if there is a risk of deposits replace the silencer at the outlet by a small flange connection or a hose nozzle and connect an outlet line, see **4.4.3 Outlet connection (EX)** on page 32.

Assemble the silencer

Assemble silencer



⇒ Unpack the silencer and screw it in the thread at the outlet of the pump.

Silencer assembled.

4.4.2 Vacuum connection (IN)



CAUTION

Flexible vacuum hoses can contract during evacuation.

Connected components that are not secured can cause injury or damage due to jerky movement (shrinkage) of the flexible vacuum hose. The vacuum hose can come loose.

- ⇒ Fix the vacuum hose to the connections.
- ⇒ Secure connected components.
- ⇒ Take the maximum shrinkage into account when sizing the flexible vacuum hose.

NOTICE

Foreign bodies in the suction line can damage the vacuum pump.


- ⇒ Prevent particles, liquids or contaminants from being aspirated or being able to flow back.

IMPORTANT!

- ⇒ Use a sufficiently stable vacuum hose that is designed for the required vacuum range.
- ⇒ Keep hose lines as short as possible.
- ⇒ Connect hose lines in a gas-tight manner to the vacuum pump.
- ⇒ Avoid kinks in the vacuum hose.

Connect the vacuum hose

→ Example
Vacuum connection
at the inlet IN

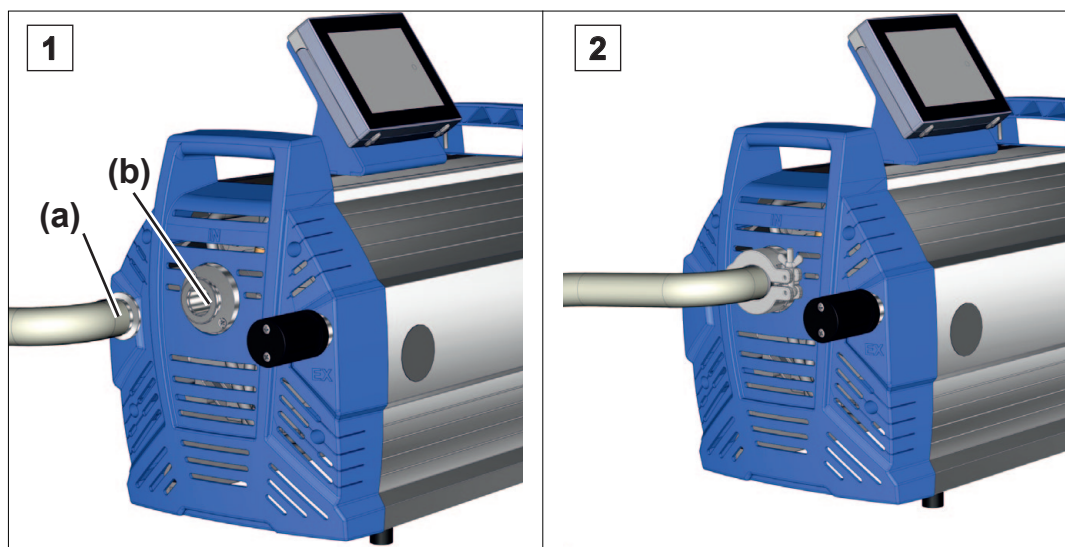

1. Take a vacuum hose (a) with a small flange connection KF DN 25.
2. Attach the vacuum hose to pump inlet (b) with a centering ring and clamping ring.




Observe the following points for optimum results:

- ⇒ Keep the vacuum line as short as you can with as large a cross-section as possible.
- ⇒ Alternatively, you can connect a vacuum hose via an adapter to the hose nozzle DN 15 mm
→ see accessories in **8.2 Ordering information on page 75**.


4.4.3 Outlet connection (EX)

As standard the diaphragm pumps are equipped with a silencer at the outlet (EX). The screwed in silencer can be replaced optionally by a small flange connection KF DN 16 or a hose nozzle DN 15 mm as outlet connections → see accessories in *8.2 Ordering information on page 75*.

	<p>WARNING</p>
	<p>Risk of bursting due to overpressure inside the outlet line.</p> <p>Inadmissibly high pressure in the outlet line can cause the vacuum pump to burst or damage seals.</p> <ul style="list-style-type: none"> ⇒ The outlet line (exhaust gas, gas outlet) must always be clear and non-pressurized. ⇒ Always route the exhaust gas hose with a fall or take measures to prevent condensate from flowing back into the vacuum pump. ⇒ Observe the maximum admissible pressures and pressure differences.

Modifying the outlet connection (optional)

→ Example
Assemble a hose nozzle at the outlet

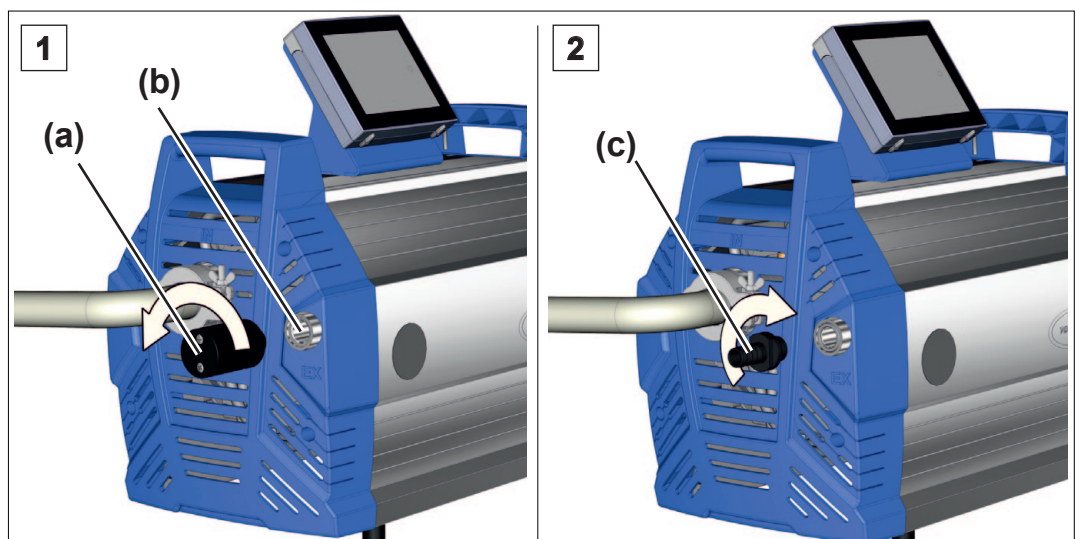


1. Unscrew the silencer (a) from the outlet connection (EX) (b).
2. Assemble a small flange connection KF DN 16 or a hose nozzle DN 15 mm (c) at the outlet connection, thread 1/2".

Connect the exhaust gas hose (optional)


Connect the exhaust gas hose to the pump, either via the small flange connection KF DN 16 or via the hose nozzle DN 15 mm. The following example describes the connection via hose nozzle.

→ Example
Exhaust gas connection at the outlet EX

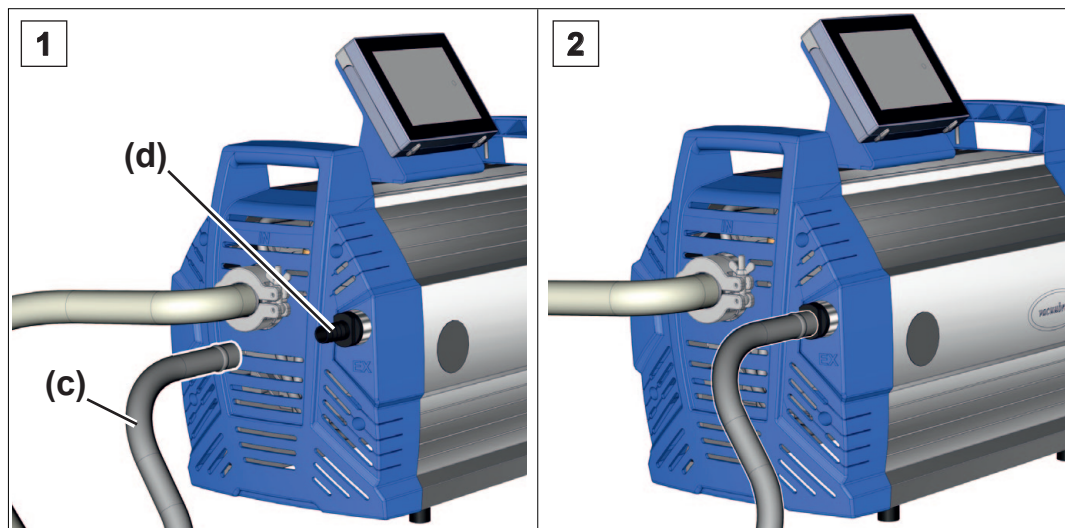

1. Take a vacuum hose (c), d_i 15 mm.
2. Push the exhaust gas hose (d) onto the hose nozzle and lay the hose, if necessary, in a fume hood. If necessary fix the outlet hose, e.g., with a hose clip.

4.4.4 Connect venting valve (optional)


	DANGER
	<p>Risk of explosion due to venting with air.</p> <p>Depending on the application, venting can cause explosive mixtures to form or other hazardous situations to arise.</p> <ul style="list-style-type: none">⇒ Never vent processes with air which could form an explosive mixture.⇒ In the case of flammable substances, use only inert gas for venting, e.g., nitrogen (max. 1.2 bar/900 Torr abs.).

There is *no* venting connection at the diaphragm pump, at the controller or at the enclosed **VACUU·SELECT® Sensor**. You can connect different venting valves, e.g., the valve **VB M-B**, via **VACUU·BUS®** directly to the controller, though. → see accessories in *8.2 Ordering information on page 75*.

4.4.5 Electrical connection

Electrical connection of the pump

→ Example
Electrical connection
diaphragm pump

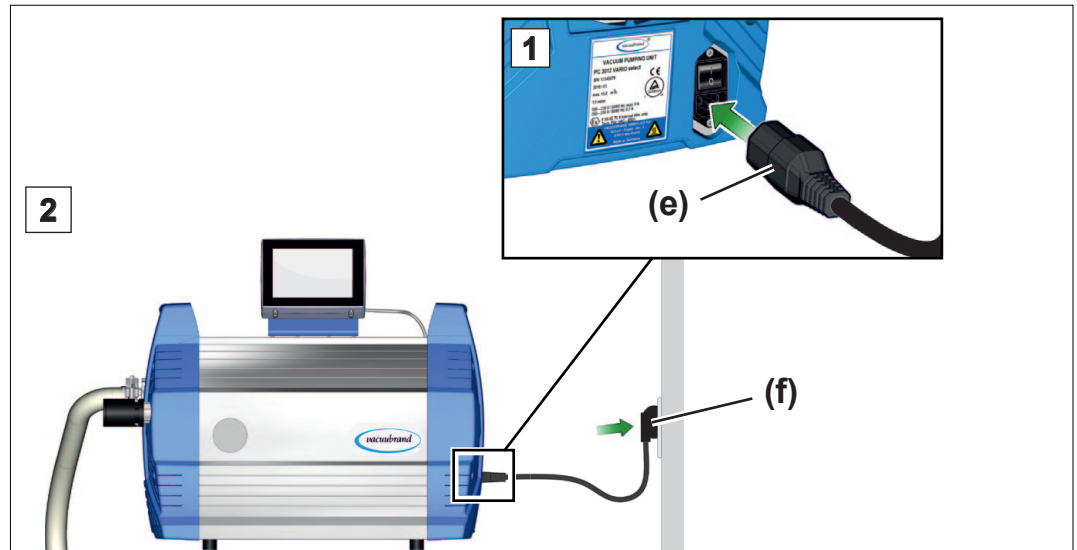

1. Plug the connector (e) on the power cable into the power connection of the vacuum pump.
2. Plug power plug (f) into the power outlet.
 - Vacuum pump electrically connected.

IMPORTANT!

⇒ Lay the power cable such that it cannot be damaged by sharp edges, chemicals, or hot surfaces.

Power connections with country code

Diagrams of
standard power
connections with
grounding contact

1		2		3		4																	
5		6		7		<table border="1"> <tr> <td>1</td> <td>IL</td> <td>5</td> <td>US</td> </tr> <tr> <td>2</td> <td>UK</td> <td>6</td> <td>CEE</td> </tr> <tr> <td>3</td> <td>CN</td> <td>7</td> <td>CH</td> </tr> <tr> <td>4</td> <td>IND</td> <td></td> <td></td> </tr> </table>		1	IL	5	US	2	UK	6	CEE	3	CN	7	CH	4	IND		
1	IL	5	US																				
2	UK	6	CEE																				
3	CN	7	CH																				
4	IND																						

The vacuum pump is delivered ready for use with the appropriate power plug.

IMPORTANT!


- ⇒ Use the power plug which fits your power supply.
 - ⇒ Do not use multiple sockets connected in series as the power connection.
-

5 Commissioning (operation)

5.1 Switch on

Switch pump on

Switch pump on

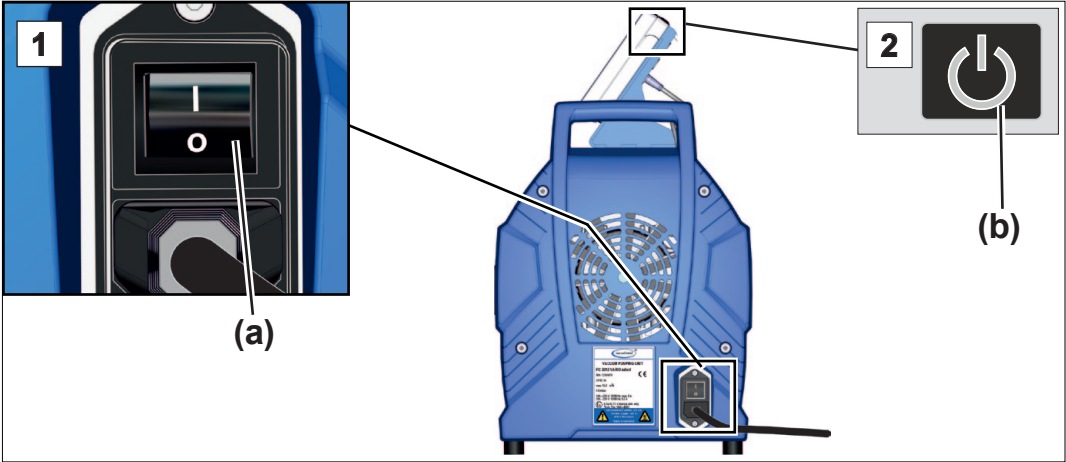


1. Switch rocker switch (a) on – switch position I.
2. Press ON/OFF button (b) on the controller.
 - The start screen is displayed.
 - After approx. 30 seconds, the process screen appears with the operating elements in the controller display.

5.2 Operation

Operation with vacuum controller


Apart from the chapters Switch on and Switch off, this manual describes the mechanical structure of the diaphragm pumps.

Operation of the vacuum controller and its functions are described in the separate **VACUU-SELECT** manual.

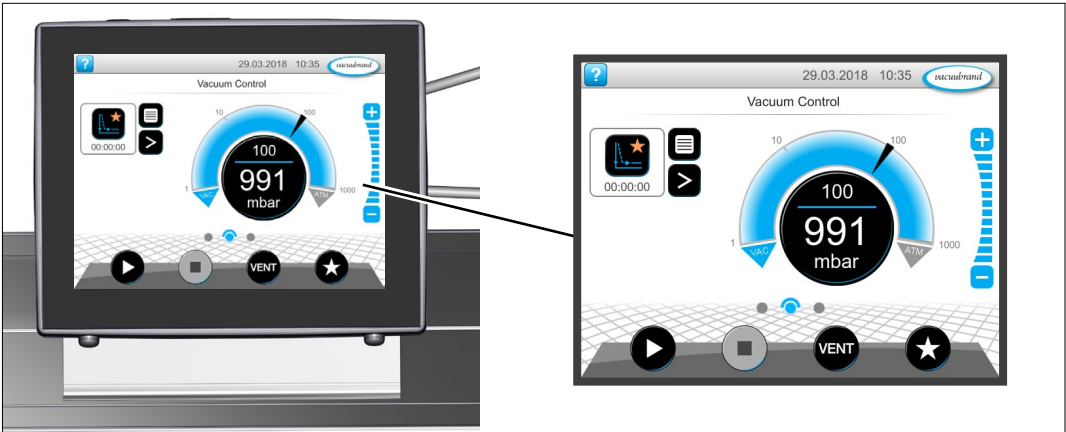

Process screen

Process screen vacuum controller

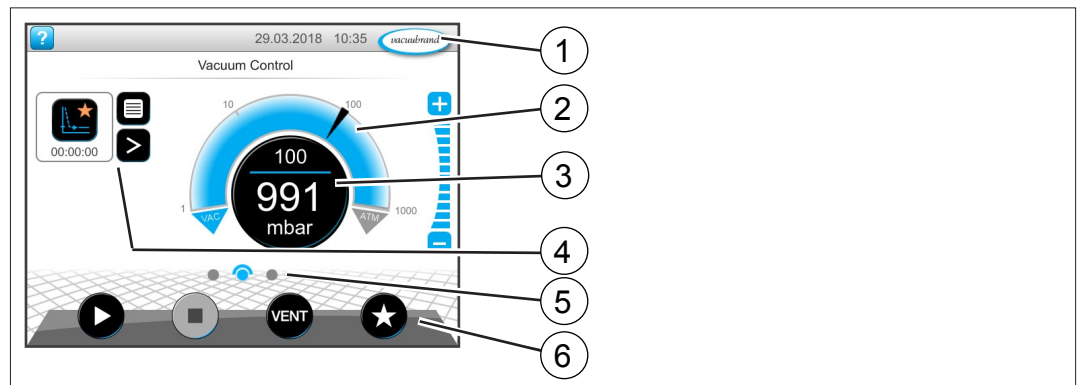

- 1 Status bar
- 2 Analogue pressure display – pressure curve
- 3 Digital pressure reading – pressure value (target value, actual value, pressure unit)
- 4 Process screen with context features
- 5 Screen navigation
- 6 Operating elements for control

Operating elements

Vacuum controller operating elements


Button		Function
active	locked	
		Start ▶ Start application – only available on the process screen.
		Stop ▶ Stop application – always possible.
		VENT – vent the system (option) ▶ Press button < 2 sec = vent briefly; control continues.
		▶ Press button > 2 sec = vent to atmospheric pressure; vacuum pump is stopped. ▶ Press button during venting = venting is stopped.
		Favorites ▶ View Favorites menu.

* Button is only displayed if venting valve is connected or activated.

5.2.1 Operation (→ see description of controller)


Start the vacuum controller

Start

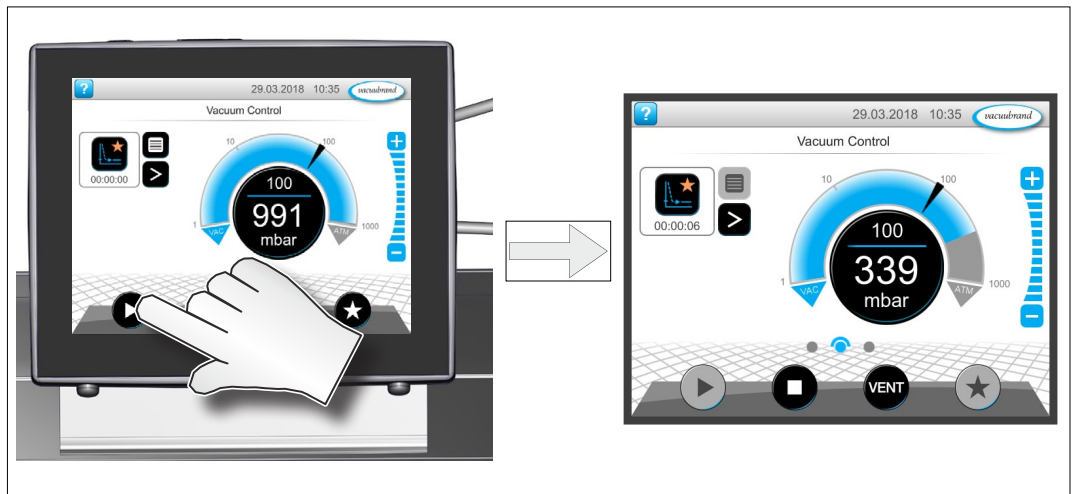


Stop the vacuum controller

Stop

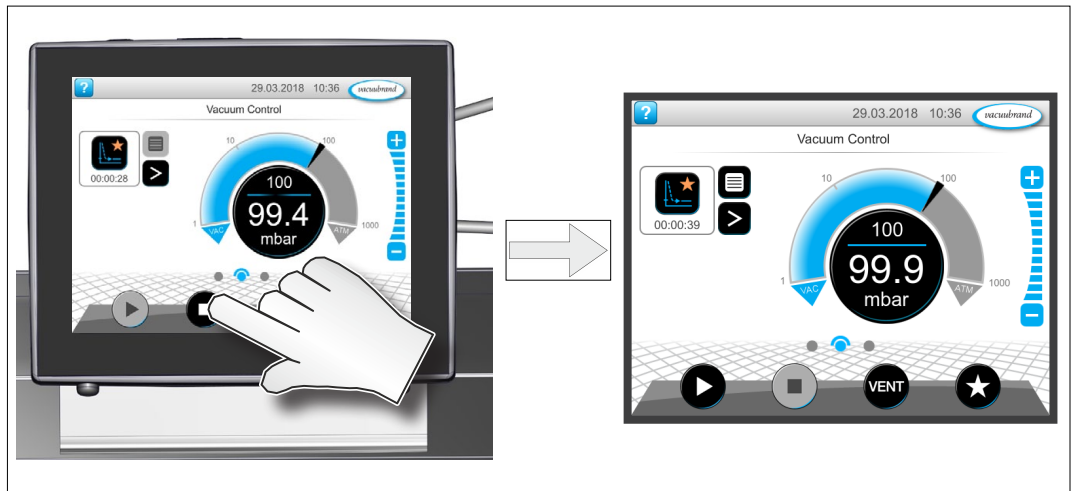

Venting (option)

Venting

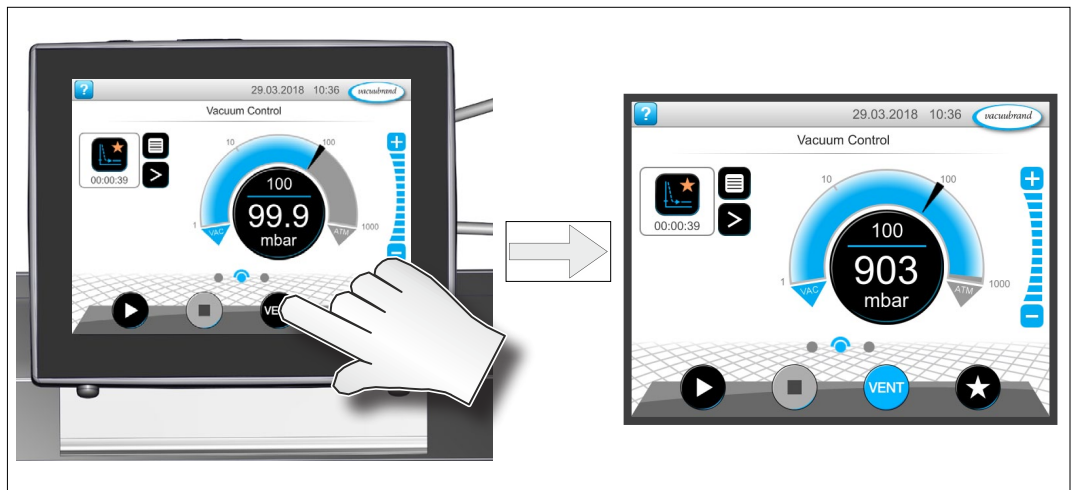


5.3 Shutdown (switch off)

Take the pump out of operation

Switch pump off

1. Stop the process.
2. Disconnect the pump from the apparatus.
3. Let the vacuum pump run on for about 30 minutes, with open inlet.
 - Condensate and media residues will be flushed out of the vacuum pump.

IMPORTANT! ⇒ Prevent deposits and rinse condensate out of the pump.

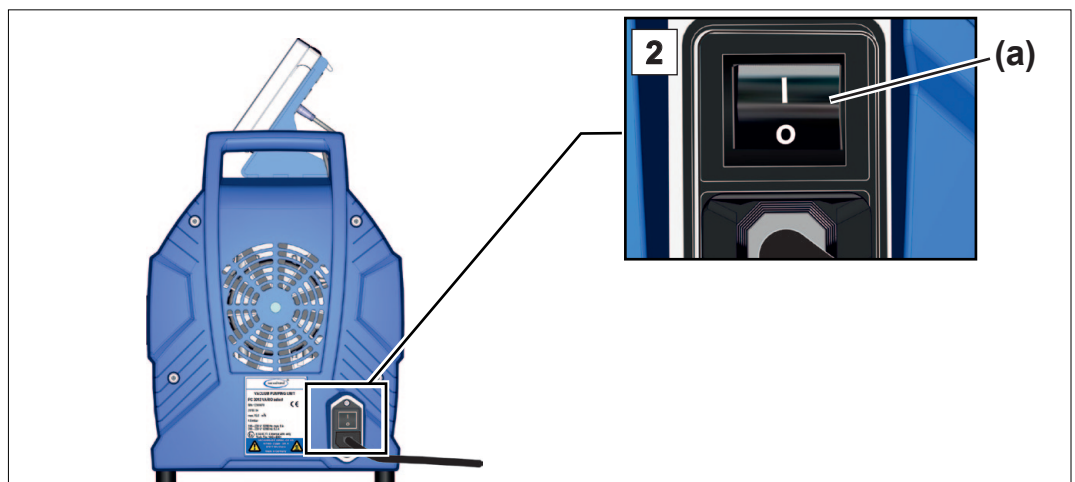

4. Switch rocker switch (a) off – switch position 0.
 - Pump switched off.
5. Check the pump for dirt and damage.

5.4 Storage

Store the vacuum pump

1. Clean the vacuum pump if dirty.
2. Recommendation: Perform a preventive maintenance before storing the vacuum pump. This is especially important if it ran more than 15,000 operating hours.
3. Close the suction and outlet lines, e.g., with the transport caps.
4. Package the vacuum pump such that it is protected from dust; enclose desiccants if necessary.
5. Store the vacuum pump in a cool, dry location.

IMPORTANT!

If damaged parts are stored for operational reasons, these should be clearly identified as **not ready for use**.


6 Troubleshooting

6.1 Technical support

⇒ To identify errors and potential remedies, please refer to the troubleshooting table *Error – Cause – Remedy*.

Technical support

For technical assistance or errors for which you require additional support, please contact your local distributor or our [Service Department](#)¹.

 Operate the machine only when it is in proper working condition.

- ⇒ Observe the recommended maintenance intervals to ensure a fully functional system.
- ⇒ Send defective devices to our Service Department or your local supplier for repair!

6.2 Error – Cause – Remedy

Error – Cause – Remedy

Error	▶ Possible cause	✓ Remedy	Personnel
Readings deviate from the reference standard	<ul style="list-style-type: none"> ▶ Vacuum sensor dirty. ▶ Moisture in the sensor. ▶ Sensor defective. ▶ Sensor measures incorrectly. 	<ul style="list-style-type: none"> ✓ Clean sensor measuring chamber. ✓ Allow sensor measuring chamber to dry, e.g., by pumping. ✓ Calibrate sensor with reference gauge. ✓ Replace defective components. 	Specialist
Sensor does not pass on measured value	<ul style="list-style-type: none"> ▶ No voltage applied. ▶ VACUU·BUS plug-in connection or cables defective or not connected. 	<ul style="list-style-type: none"> ✓ Check VACUU·BUS plug-in connection and cables to the controller. 	Operator
	<ul style="list-style-type: none"> ▶ Sensor defective. 	<ul style="list-style-type: none"> ✓ Replace defective components. 	Specialist

¹ -> Phone: +49 9342 808-5660, fax: +49 9342 808-5555, service@vacuubrand.com

Error – Cause –
Remedy




Error	▶ Possible cause	✓ Remedy	Personnel
Venting valve (optional) does not switch	<ul style="list-style-type: none"> ▶ No voltage applied. ▶ VACUU·BUS plug-in connection or cables defective or not connected. ▶ Venting valve dirty. 	<ul style="list-style-type: none"> ✓ Check VACUU·BUS plug-in connection and cables to the controller. ✓ Clean venting valve. ✓ Perform component detection in VACUU·SELECT – see: Main menu/ Administration/ VACUU·BUS. ✓ If necessary, use another external venting valve. 	Specialist
Vacuum pump does not start	<ul style="list-style-type: none"> ▶ Overpressure in the outlet line. ▶ Condensation in the vacuum pump. 	<ul style="list-style-type: none"> ✓ Open the outlet line, check the silencer. ✓ Ensure clear passage. 	Operator
	<ul style="list-style-type: none"> ▶ Pump switched off. ▶ Power plug not correctly plugged in or pulled out. ▶ VACUU·BUS plug-in connection or cables defective or not connected. 	<ul style="list-style-type: none"> ✓ Switch pump on using rocker switch. ✓ Check power supply and cable. ✓ Check VACUU·BUS plug-in connection and cables to the controller. 	Operator
	<ul style="list-style-type: none"> ▶ Motor overloaded. ▶ Thermal protection has been triggered. 	<ul style="list-style-type: none"> ✓ Allow the motor to cool down. ✓ Clear error manually: <ul style="list-style-type: none"> → Unplug pump from the power supply → Eliminate cause of error → Switch pump back on 	Specialist

Error – Cause – Remedy

Error	▶ Possible cause	✓ Remedy	Personnel
No or very little suction power	▶ Leak in the suction line or in the apparatus.	✓ Check suction line and apparatus for leaks.	Operator
	▶ Vacuum line too long or cross-section too small.	✓ Use a shorter vacuum line with a larger cross-section.	Operator
	▶ Condensate inside the vacuum pump.	✓ Allow vacuum pump to run for a few minutes with the suction nozzle open.	Operator
	▶ Deposits inside the vacuum pump.	✓ Clean and check pump heads.	Specialist
	▶ Diaphragms or valves defective.	✓ Replace diaphragms and valves.	Specialist
	▶ High level of vapor generated in the process.	✓ Check process parameter.	Specialist
No display	▶ Pump switched off. ▶ Power plug not correctly plugged in or pulled out. ▶ VACUU·BUS plug-in connection or cables defective or not connected. ▶ Controller switched off or defective.	✓ Switch pump on using rocker switch. ✓ Switch on controller. ✓ Check power supply and cable. ✓ Check VACUU·BUS plug-in connection and cables to the controller.	Operator
		✓ Replace defective components	Specialist

Error	▶ Possible cause	✓ Remedy	Personnel
Loud operating noises	<ul style="list-style-type: none">▶ No silencer or hose mounted at the outlet.▶ Outlet line open.	<ul style="list-style-type: none">✓ Check silencer or hose and install correctly.✓ Check outlet line connections.✓ Connect the outlet line to an extraction system or fume hood.	Operator
	<ul style="list-style-type: none">▶ Ball bearing defective.	<ul style="list-style-type: none">✓ Service the vacuum pump and replace defective parts or send in the device.	Specialist

7 Cleaning and maintenance

 	<p>WARNING</p> <p>Danger due to electrical voltage.</p> <ul style="list-style-type: none"> ⇒ Switch the device off before cleaning or maintenance work. ⇒ Unplug the power plug from the socket.
	<p>Risk from contaminated parts.</p> <p>Pumping hazardous media can result in hazardous substances adhering to internal parts of the pump.</p> <ul style="list-style-type: none"> ⇒ Wear your personal protective equipment, e.g., protective gloves, eye protection and, if necessary, respiratory protection. ⇒ Decontaminate the vacuum pump before opening it. If necessary have decontamination carried out by an external service provider. ⇒ Take safety precautions according to your instructions for handling hazardous substances.

<p>NOTICE</p>	
<p>Damage possible if work is performed incorrectly.</p> <ul style="list-style-type: none"> ⇒ Have maintenance work performed by a trained professional or at least by a trained person. ⇒ Recommendation: Before carrying out maintenance for the first time, please read through all the instructions to get an overview of the required service work. 	

7.1 Information on service work

Recommended maintenance intervals

Maintenance intervals*	if required	15000 h
Replace diaphragms		X
Replace valves		X
Replace O-rings		X
Clean or replace molded PTFE-hose	X	
Clean the vacuum pump	X	

* Recommended maintenance interval after hours of operation and under normal operating conditions; depending on the environment and area of application, we advise performing cleaning and maintenance as needed.

Recommended aids

→ Example
Recommended aids
for cleaning and
maintenance



Protective gloves

IMPORTANT!

⇒ Always wear your personal protective equipment when performing activities which may bring you into contact with hazardous substances.

Tools needed for maintenance

→ Example
Tool



Nr	Tool	Size
1	Service kit	
	Service kit MD 12 / MV 10 VARIO select #20696827	1x
	<i>or</i>	
	Service kit ME 16 VARIO select #20696819	1x
2	Diaphragm wrench #20636554	SW66
3	Flat nose pliers	
	Close hose clamps	
4	Flat-head screwdriver	
	Open hose clamps	Size 1
5	Phillips screwdriver	
	Screw fittings, controller base	Size 1
	Screw fittings, distributor, outlet holder	Size 2
6	Hex key	
	Screw fittings, side panels	Size 5
	Screw fittings, housing cover	Size 5
	Screw fittings, housing sections with handle	Size 4
	Screw fittings, side panel retaining plates	Size 4
7	Torque wrench, adjustable 2–10 Nm	

7.2 Cleaning

IMPORTANT!

This chapter does not contain descriptions for decontamination of the product. This chapter describes simple measures for cleaning and care.

⇒ Before cleaning, switch off the diaphragm pump.

7.2.1 Diaphragm pump

Clean the surfaces



Clean dirty surfaces with a clean, slightly damp cloth. We recommend using water or mild soapy water to moisten the cloth.

7.2.2 Clean or replace molded PTFE hoses

Maintenance provides the opportunity to check the components of the diaphragm pump, including the hoses.

⇒ Clean the inside of very dirty molded hoses, e.g., using a pipe cleaner or similar.

⇒ Replace brittle and defective molded hoses.

7.2.3 Clean or replace the controller

During maintenance, the controller can be disconnected and removed.

Clean the surfaces



⇒ Clean dirty surfaces with a clean, slightly damp cloth. We recommend using water or mild soapy water to moisten the cloth.


⇒ Reattach the controller after maintenance work has been completed.

7.3 Diaphragm pump maintenance

7.3.1 Maintenance items

Servicing positions

→ Example
Diaphragm pump,
front,
semi-transparent
view

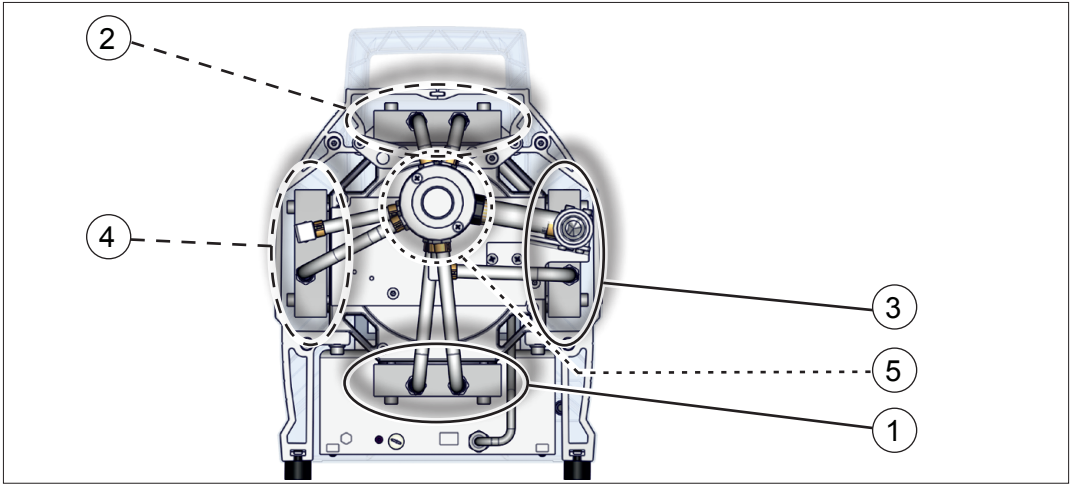

Meaning

Maintenance items and sequence	
1	Bottom pump head pair
2	Top pump head pair
3	Right pump head pair
4	Left pump head pair
5	Suction/pressure distributor (only MV 10 / MD 12)

i Straightforward maintenance due to split work steps. Observe the recommended sequence of maintenance steps according to the table:

- ⇒ On one pump head pair, first replace the diaphragms.
- ⇒ Then change the inlet/outlet valves.
- ⇒ Repeat these steps on the next pump head pair.
- ⇒ Then replace the O-ring and the pressure relief valve in the suction/pressure distributor.


7.3.2 Preparation

Disassemble the controller and base


→ see also chapter: *4.3 Controller base on page 28*

Disassemble the device and housing sections

Disassemble the front housing section

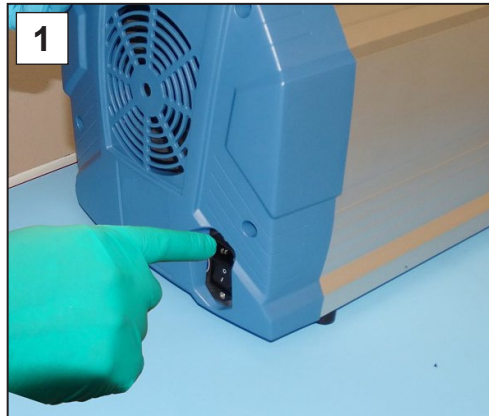


1. Switch the diaphragm pump off and unplug the power plug.

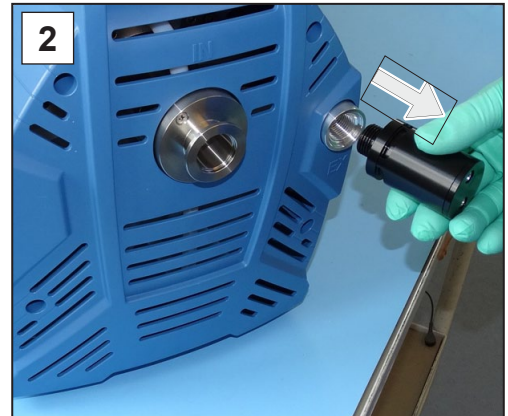


2. Unscrew the silencer from the outlet.




3. Unscrew the 4 screws from the front housing section; hex key size 4.

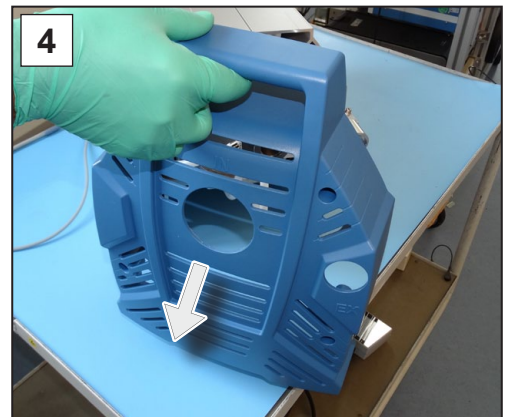

4. Remove the housing section and set it aside.

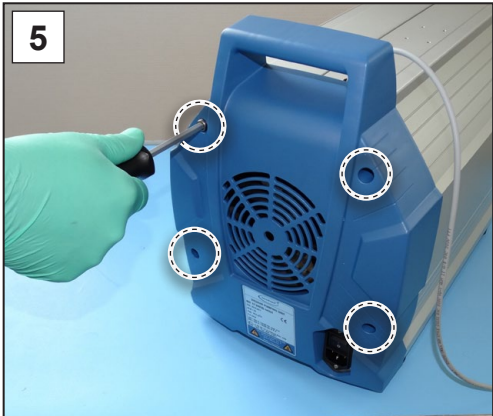


5. Unscrew the 4 screws from the rear housing section; hex key size 4.

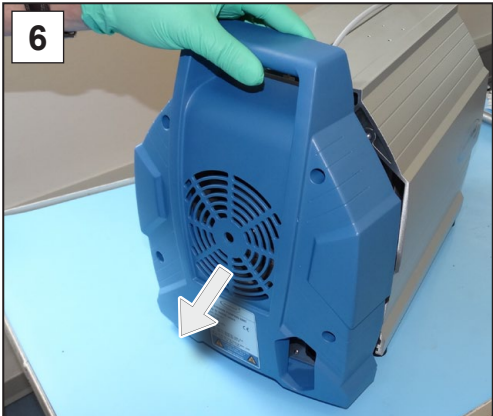


6. Remove the housing section and set it aside.

Remove the side panel

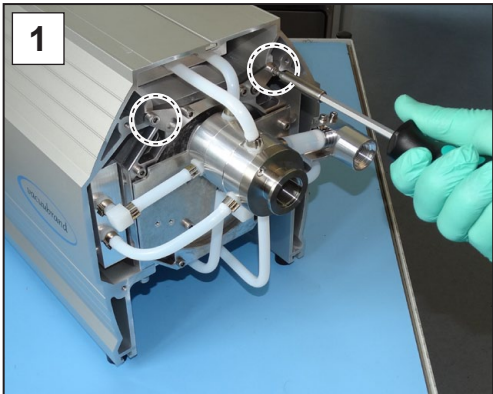


1. Unscrew the 2 outer screws from the side panel retaining plate at the front; hex key size 4

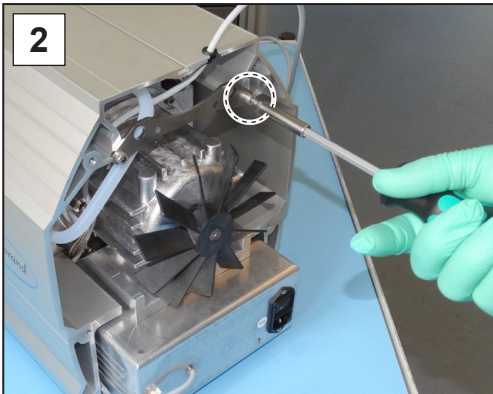


2. Unscrew the right screw from the side panel retaining plate at the rear; hex key size 4


Remove the left side panel

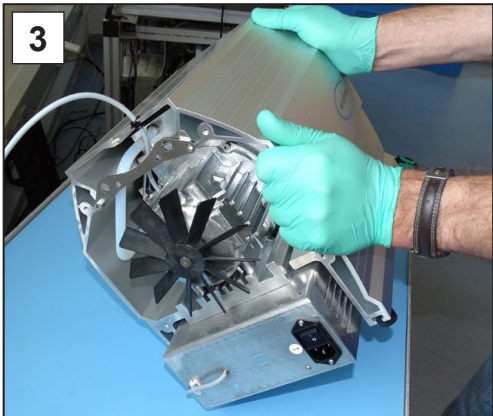


3. Place the pump carefully on its side.

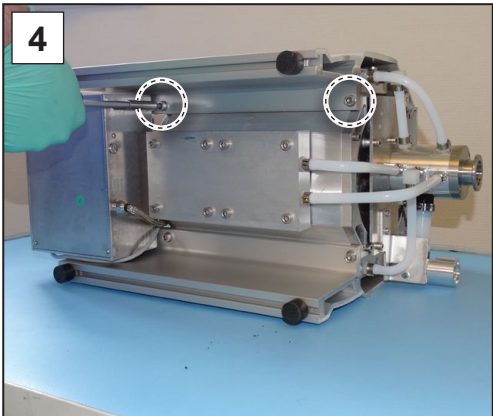

4. Unscrew the screw fittings from the side panel; hex key size 5.

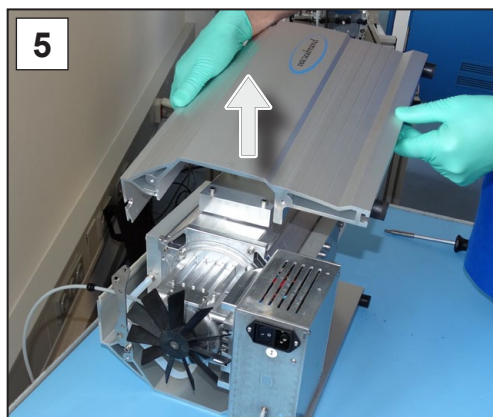


5. Lift the left side panel off the pump. Place the pump carefully on its other side.

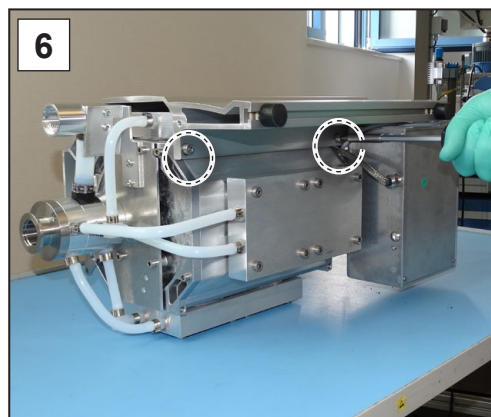

6. Unscrew the screw fittings from the side panel; hex key size 5.

Remove the right side panel

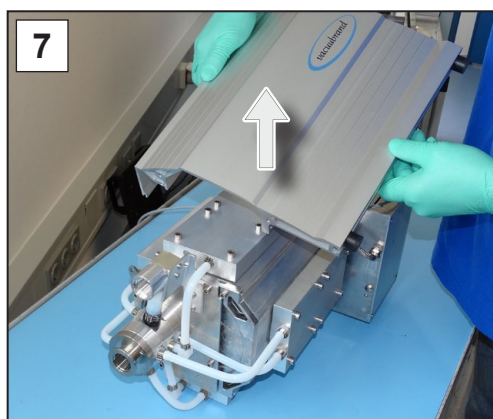


7. Lift the right side panel off the pump.

7.3.3 Replacing the diaphragms and valves

Exploded drawing of pump head (example)

Exploded-drawing
pump head

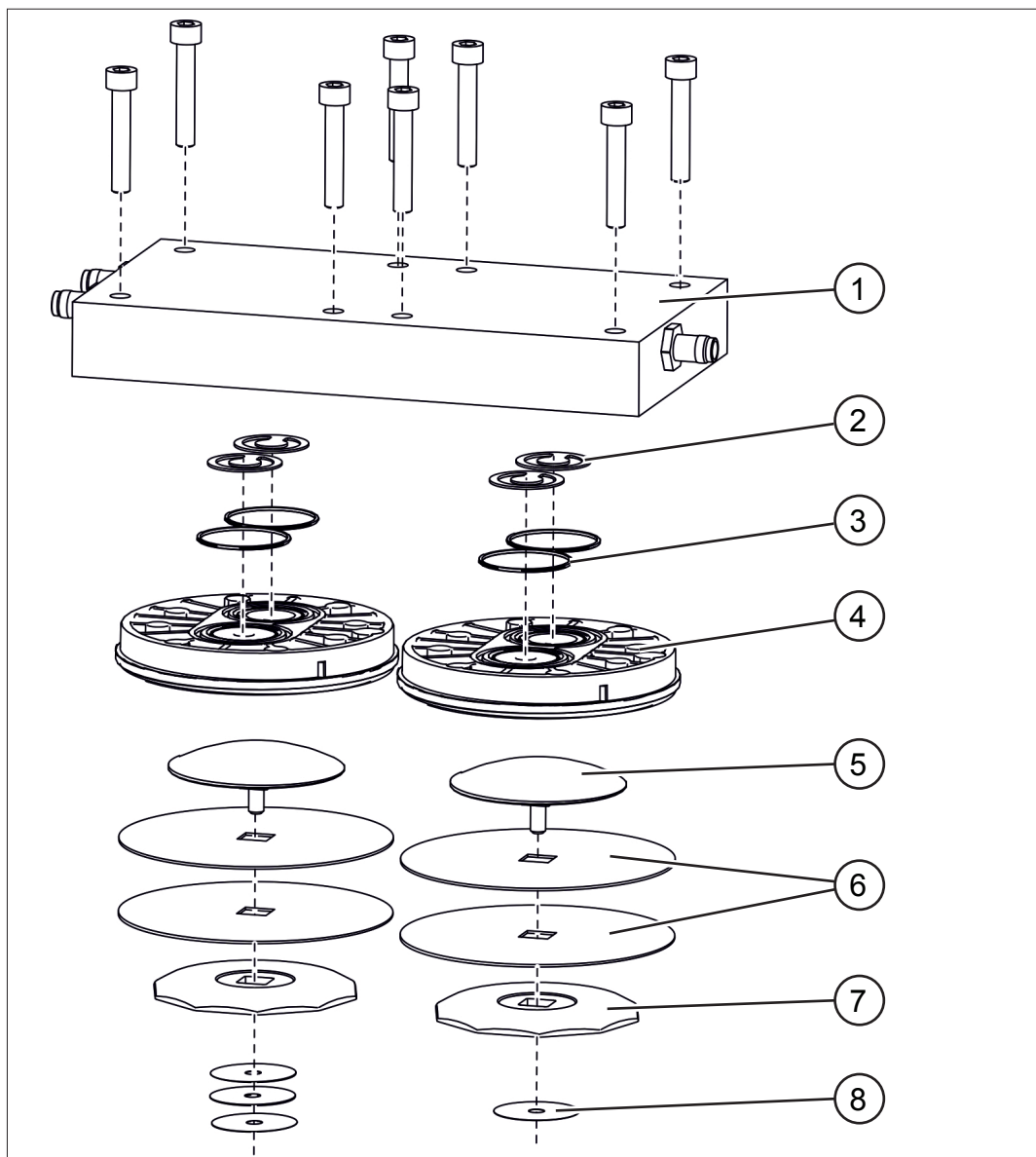

Valve maintenance


- 1 Housing cover
- 2 Valves
- 3 O-rings size 28 x 1,5

Diaphragm maintenance


- 4 Head cover
- 5 Diaphragm clamping disc with square-head screw
- 6 Double diaphragm, 2 diaphragms per pump head
- 7 Diaphragm support disc
- 8 Spacer discs, max. 4 per pump head

Bottom pump head pair


→ Example
Bottom pump head
pair

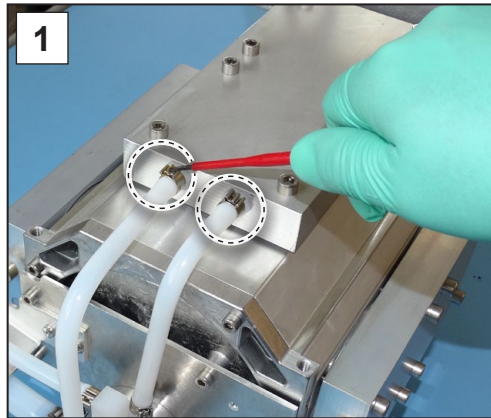


1. Turn the pump to bring the bottom pump head pair to the top. Open the hose clips on the hoses. Flat-head screwdriver size 1.

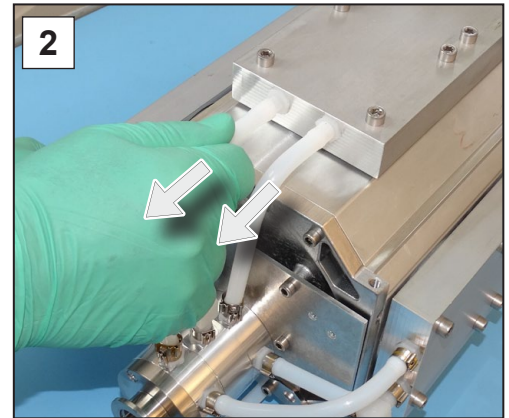


2. Pull off the molded hoses.

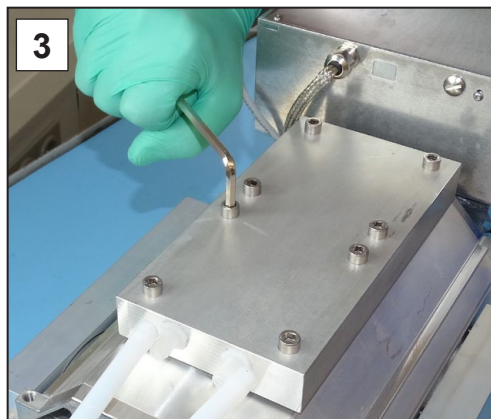

3. Unscrew the socket head screws from the housing cover. Hex key size 5.

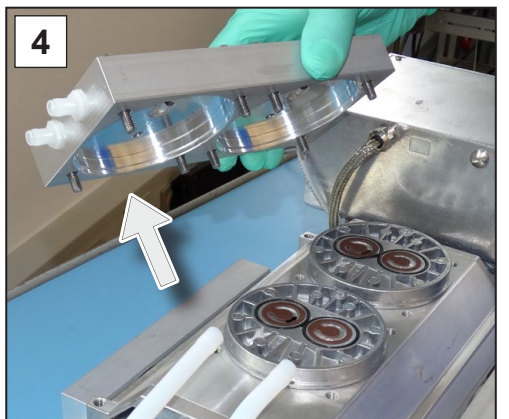


4. Remove the housing cover.

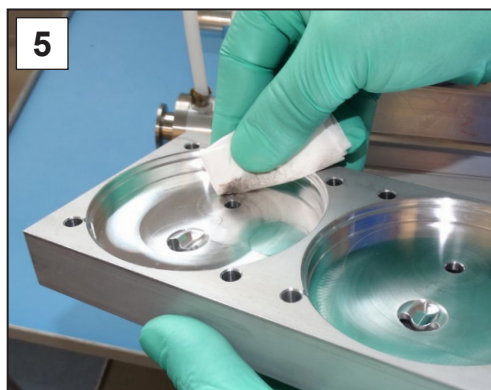

5. Check the surfaces for dirt. Clean dirty surfaces carefully.


Remove valves and O-rings

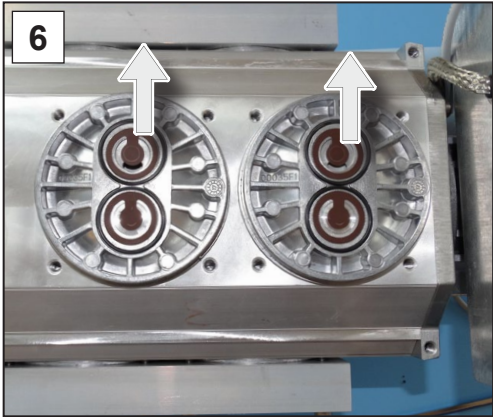

6. Carefully remove the used valves.




7. Carefully remove the used O-rings.

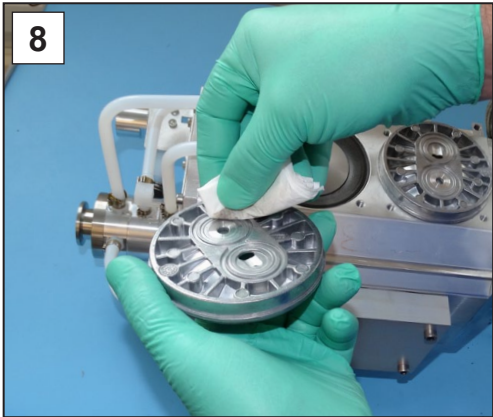

8. Remove the head covers. Check the surfaces for dirt. Clean dirty surfaces carefully.

Replace the diaphragms


→ Example
Replacing the
diaphragms

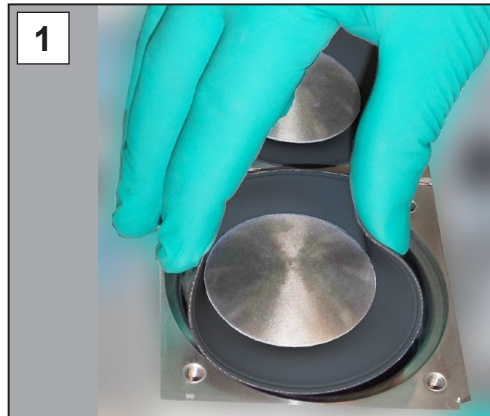

1. Lift the diaphragm upwards on either side.



2. Carefully position the diaphragm wrench on the diaphragm support disc and unscrew the assembly with the diaphragm wrench attached.

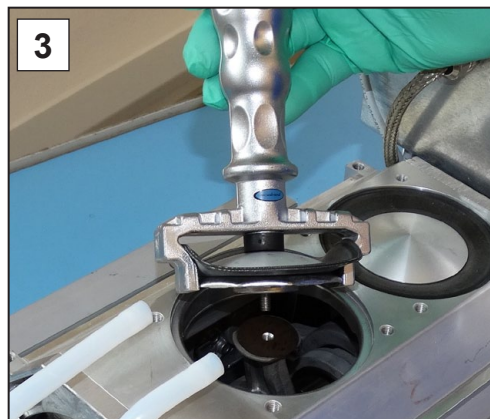

3. Lift the diaphragm, along with all the parts, out of the vacuum pump.



4. If the spacer discs adhere to the connecting rod, remove them carefully.

IMPORTANT!

- ⇒ Never drop spacer discs into the aluminum housing.
- ⇒ Keep the spacer discs. It is essential to reinsert the same number of spacer discs.


→ Example
Replacing the
diaphragms



5. Pull out the diaphragm clamping disc and remove the used diaphragm.



6. Place the new diaphragm over the square head of the clamping disc.

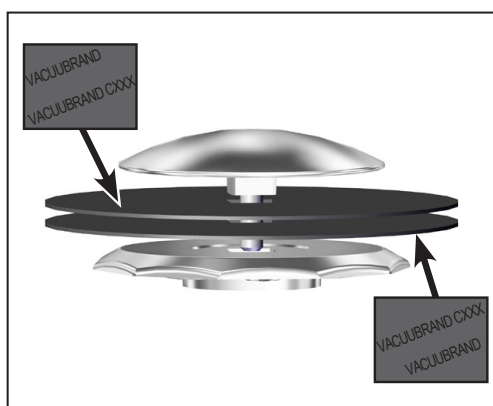

IMPORTANT!

⇒ Double diaphragm comprising 2 single diaphragms. Use the diaphragms only in pairs. The printed surfaces of the diaphragms have to face outwards.


⇒ Pay special attention to correct positioning on the square head.




7. Place all spacer discs on the thread pin.



8. Secure the diaphragm assembly inside the diaphragm wrench.

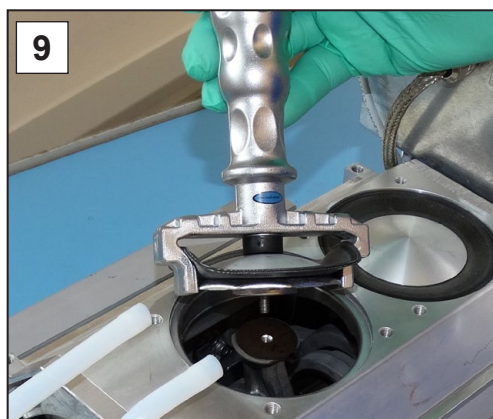


9. Hold the spacer discs firmly and place all the components carefully on the connecting rod thread.

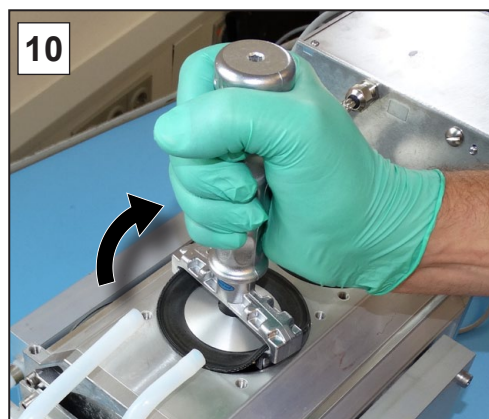

10. Initially tighten the assembly with the diaphragm wrench by hand.




11. Then position a torque wrench with socket head bit on the diaphragm wrench and tighten the assembly to 6 Nm.

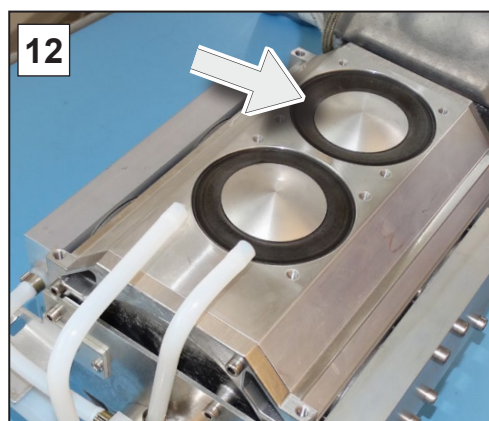

12. Repeat steps 1-11 for replacing the next diaphragm.

Insert valves


→ Example
Insert valves and
O-rings

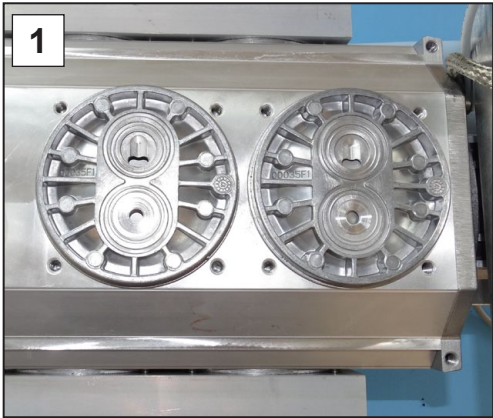


1. Place the head covers onto the diaphragms. Pay attention to the correct orientation of the head covers.




2. Insert the new O-rings into the grooves. Insert the new valves. Pay attention to the correct orientation of the valves.

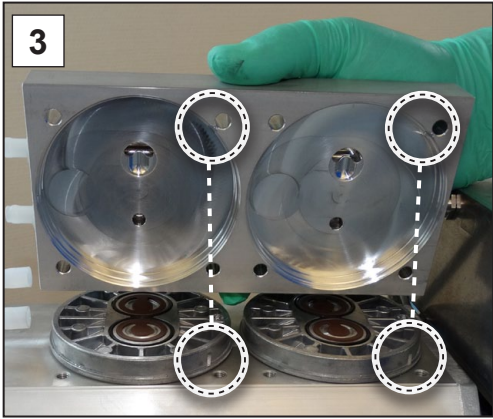


3. Take the housing cover and align the head covers. The nibs at the head covers have to lock into the notches of the housing cover.

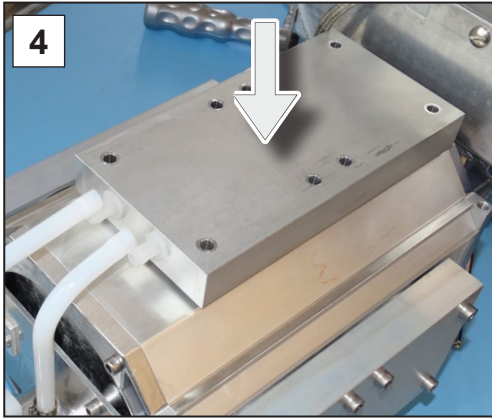


4. Put the housing cover on properly.

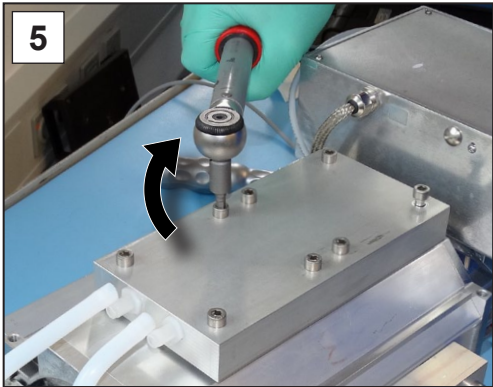


5. Screw in the screw fittings crosswise. Tighten the screw fittings with a torque wrench to 6 Nm; hex key size 5.

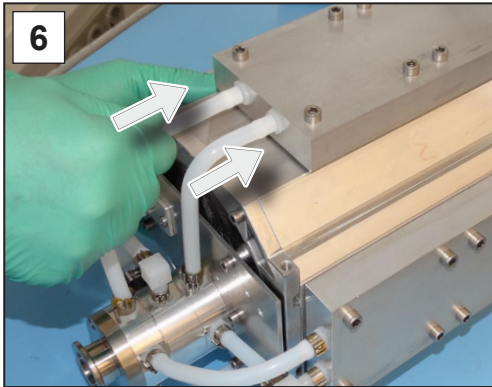

6. Slide the molded hoses back onto the hose nozzles.

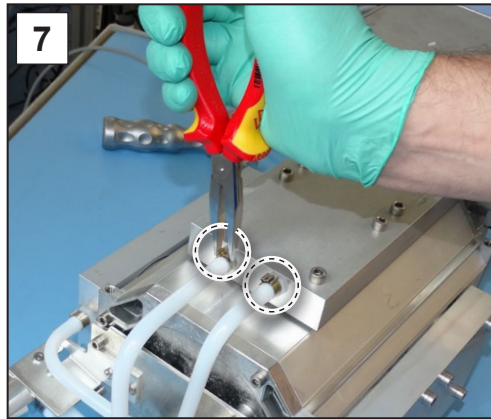

7. Secure the hose clips on the hose nozzles, e.g., with flat nose pliers.

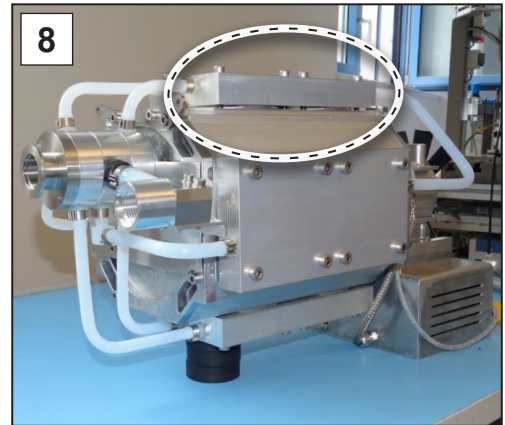

8. Turn the pump to bring the top pump head pair to the top. Support the pump, e. g., with rigid foam.


Maintain the top, right and left pump head pairs

⇒ Follow the same procedure to replace the diaphragms and valves of the top pump head pair as described for the **Bottom pump head pair, on pages 56 to 62.**


⇒ Subsequently replace one by one the diaphragms and valves of the right and left pump head pair.

Assemble the device and housing sections


Mount the side panel

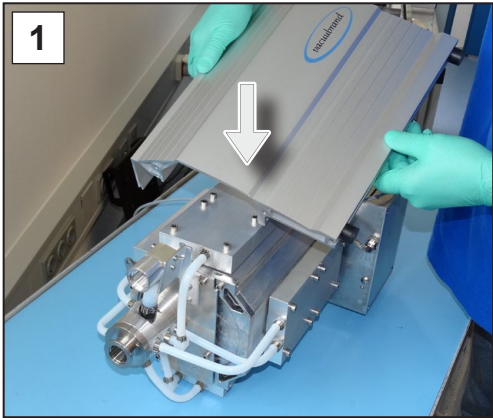


1. Place the side panel on the pump.

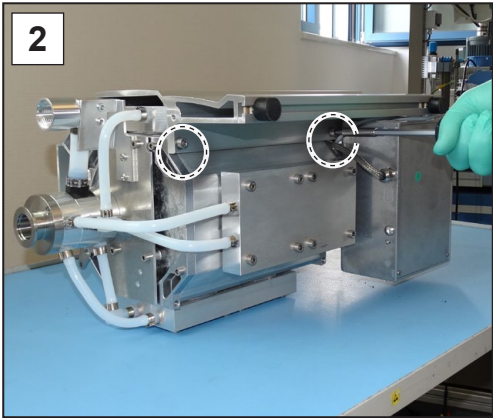


2. Wind the screw fittings into the side panel; hex key size 5.

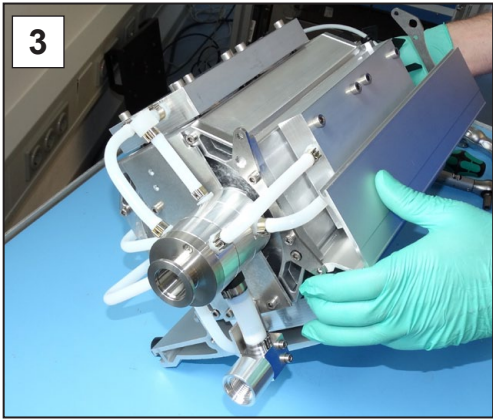


3. Turn the pump to the top.

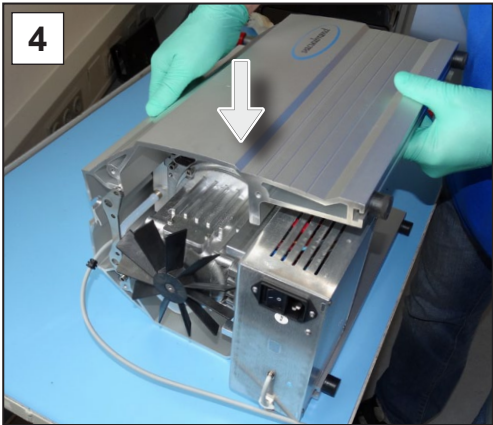


4. Place the side panel on the pump.

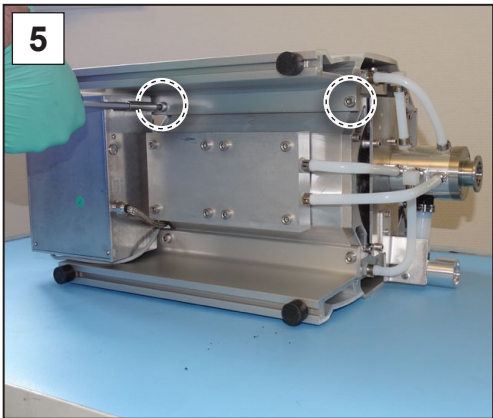


5. Wind the screw fittings into the side panel; hex key size 5.

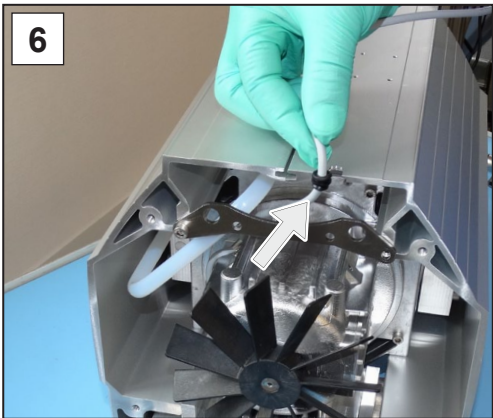


6. Secure the cable in the rear recess.

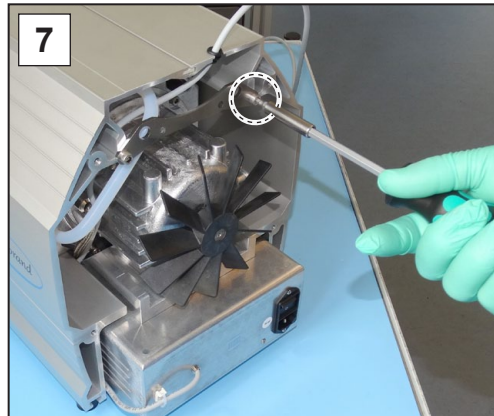

7. Wind in the screw of the side panel retaining plate at the rear; hex key size 4.

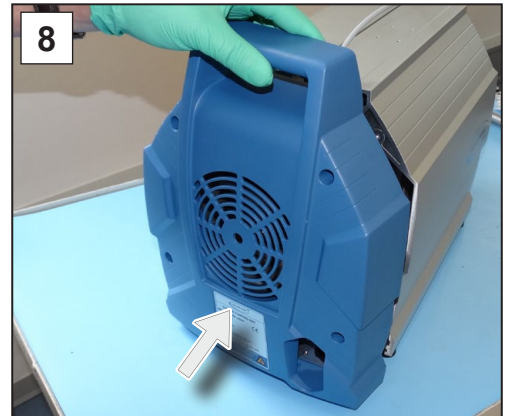

8. Place the rear housing section.

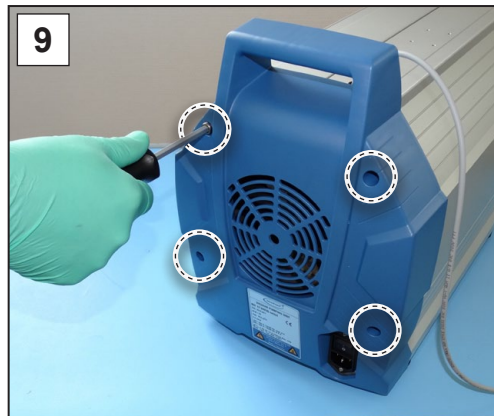

9. Wind in the screws of the housing section; hex key size 4.

⇒ Prior to mounting the front housing section the suction/pressure distributor of the pumps MD 12 and MV 10 VARIO select has to be maintained.


Suction/pressure distributor maintenance

Maintenance of suction/pressure distributor

This description only applies to diaphragm pumps MD 12 and MV 10 VARIO select.

Exploded drawing of suction/pressure distributor (example)

→ Example Pressure relief valve

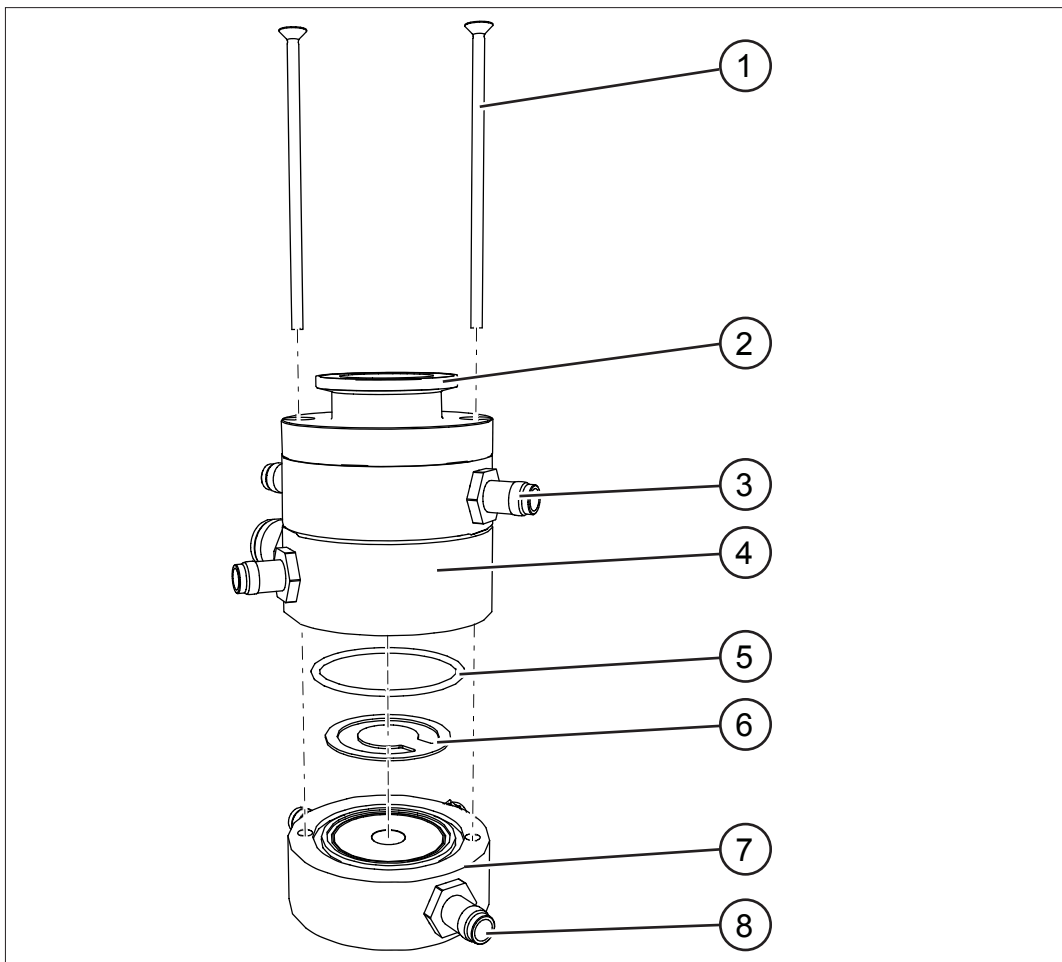


Maintenance overpressure relief valve + O-ring


1	Countersunk screw M4x80
2	Connection DN 25
3	Hose nozzle
4	Suction distributor
5	O-ring 40 x 2
6	Pressure relief valve D37
7	Pressure distributor
8	Hose nozzle

Replace pressure relief valve + O-ring


Replace pressure relief valve and O-ring

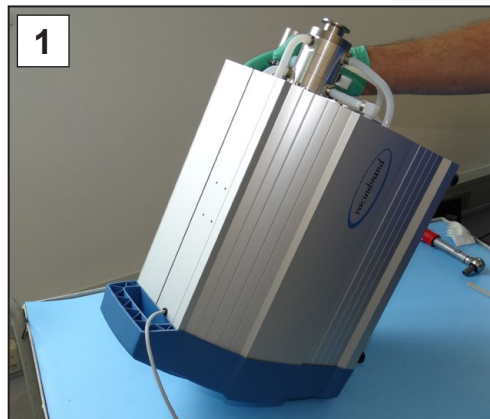


1. Place the vacuum pump on a clean, stable surface as shown.

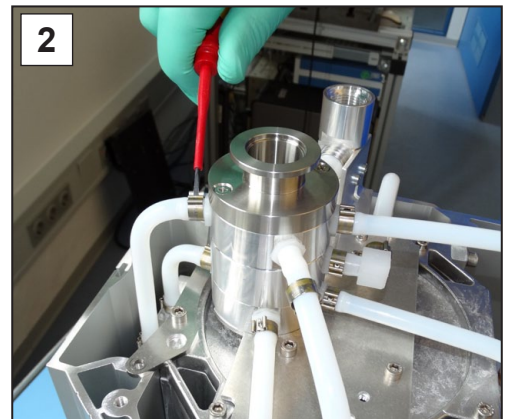


2. Only open the hose clips above the pressure distributor; flat-head screwdriver size 1.

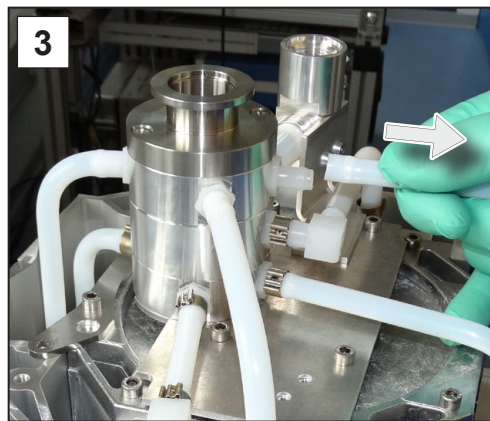


3. Remove the molded hoses one by one from the hose nozzles.

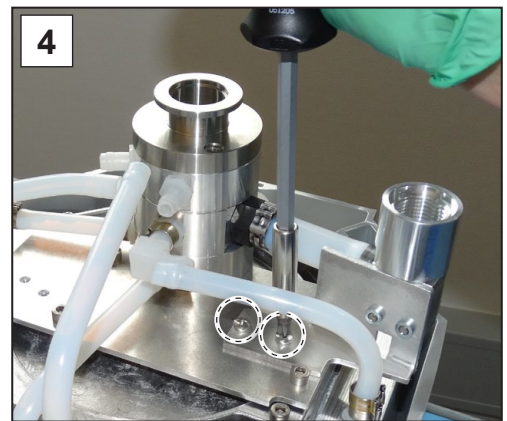


4. Unscrew the screws of the outlet holder. Phillips screwdriver size 2.

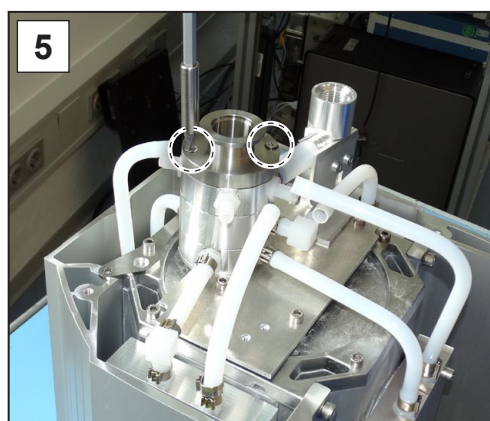


5. Unscrew the screws of the distributor. Phillips screwdriver size 2.

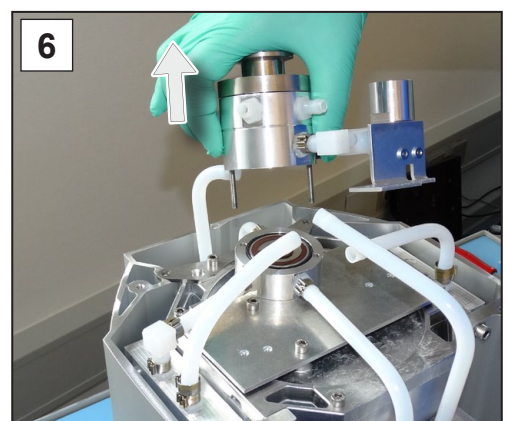


6. Remove the suction distributor together with the screws and the outlet holder and put it aside.

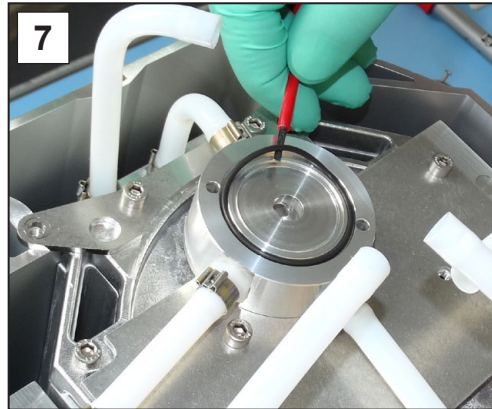


7. Carefully remove the used pressure relief valve and the O-ring, e. g., with a narrow flat-head screwdriver.

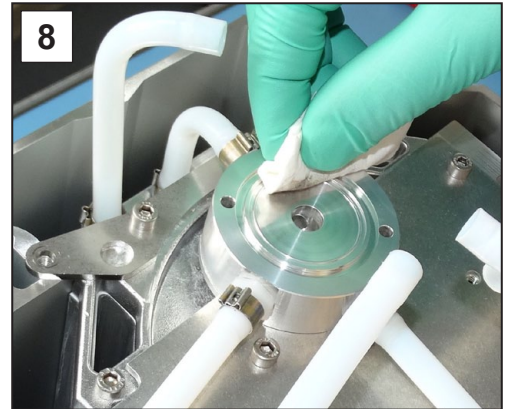


8. Clean the pressure distributor if necessary.

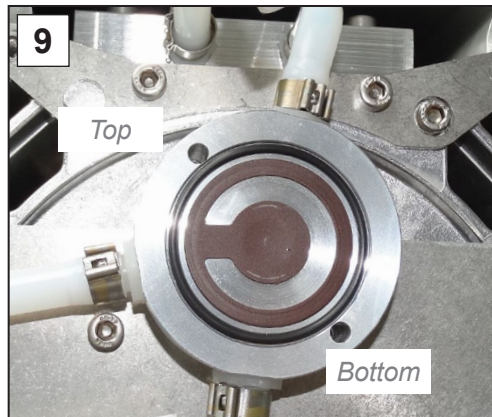

9. Place the new pressure relief valve on the clean surface. Ensure the pressure relief valve is positioned correctly on the pressure distributor. Insert the O-ring.

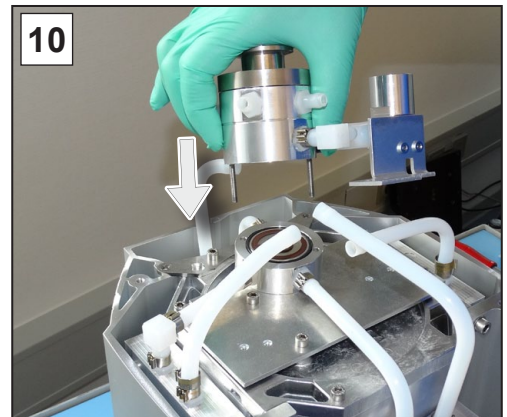


10. Place the suction distributor with screws and outlet holder onto the pressure distributor.

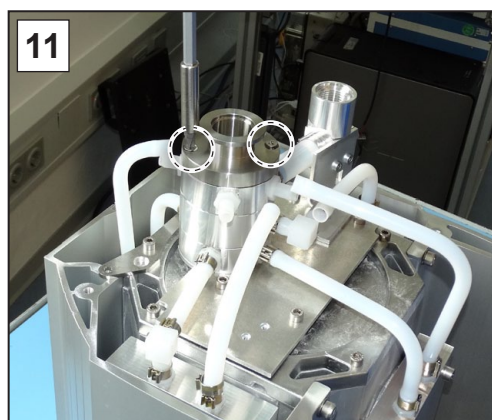


11. Wind in the screws at the distributor hand-tight; Phillips screwdriver size 2.



12. Wind in the screws at the outlet holder hand-tight; Phillips screwdriver size 2.

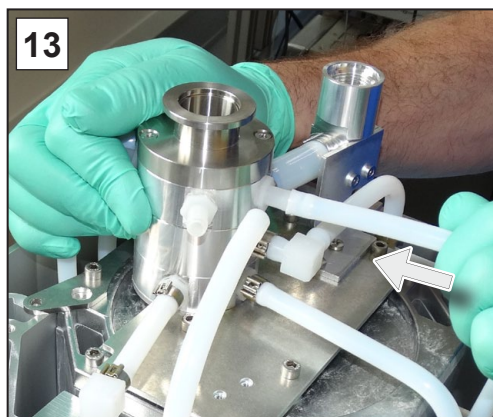

13. Push the molded hoses back into place on the hose nozzles.

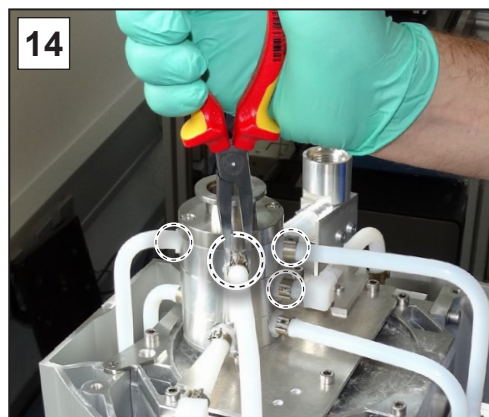


14. Secure the hose clips on the hose nozzles, e.g., with flat nose pliers.

Assemble the device and housing sections


Assemble the device and housing sections




1. Insert the bar into the groove between the side panels.

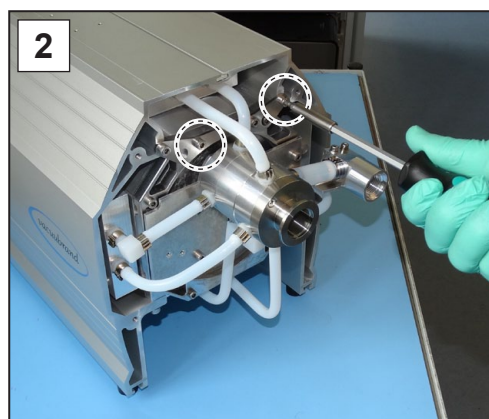


2. Wind in the 2 outer screws of the side panel retaining plate at the front; hex key size 4.

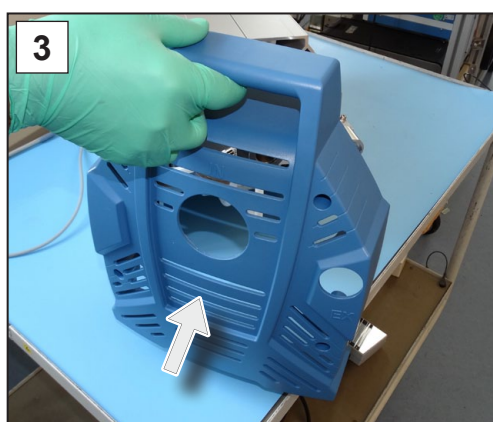

3. Place the front housing section.

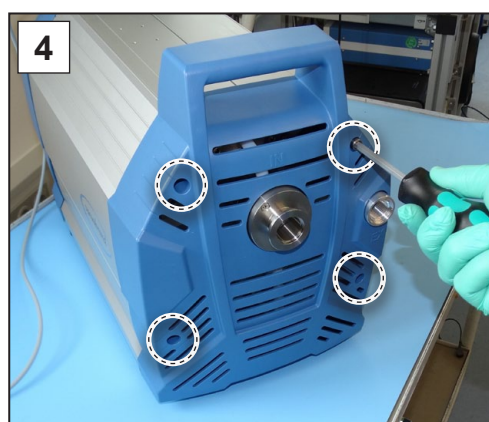

4. Wind in the 4 screws of the front housing section; hex key size 4.

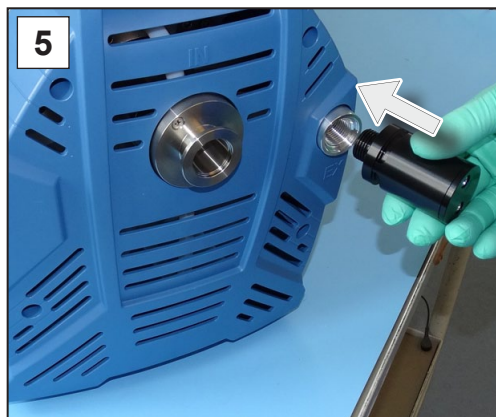

5. Screw the silencer in the thread at the outlet..



6. Secure the controller on the diaphragm pump and connect all cables.



7. Plug in the power plug.

If maintenance work has been completed in full:

- ⇒ Connect the hoses for operation.
- ⇒ Connect the diaphragm pump to the power supply.
 - Diaphragm pump is ready to be returned to use.

If not reconnected:

- Diaphragm pump is ready for storage.

8 Appendix


8.1 Technical information

Diaphragm pump series	
ME 16 VARIO select	MD 12 VARIO select
MV 10 VARIO select	

8.1.1 Technical data

Technical data

Ambient conditions		(US)
Ambient temperature, max.	10–40 °C	50–104°F
Working temperature	10–40 °C	50–104°F
Storage/transport temperature	-10–60 °C	14–140°F
Max. altitude	2000 m above NHN	6562 ft above sea level
Relative humidity	30–85 %, non condensing	
Protection class	IP 40 / IK 08	

Operating conditions		(US)
Maximum admissible media temperature (gas), non-explosive atmosphere:		
Short term	80 °C	176°F
Continuous operation	40 °C	104°F
ATEX conformity	II 3/- G IIC T3 X internal atm. only	
Maximum admissible media temperature (gas)  atmosphere:		
Short term	40 °C	104°F
Continuous operation	40 °C	104°F

Connections	
Vacuum, inlet	Small flange KF DN 25
Exhaust gas, outlet EX	Silencer / thread 1/2"
Cold device plug	+ power supply CEE, CH, CN, UK, IN, US
Plug-in connector	VACUU·BUS®

Technical data

Electrical data		(US)
Nominal voltage	200-230 VAC	100-120 VAC
Nominal frequency	50 Hz / 60 Hz	50 Hz / 60 Hz
Nominal current	6,3 A	8 A
Power, max.	1 kW	
Interface	VACUU·BUS®	
Power cable	2 m	
Vacuum data		(US)
ME 16 VARIO select		
Max. pumping speed	20 m ³ /h	11.8 cfm
Ultimate vacuum, abs.	70 mbar	53 Torr
Number of cylinders/stages	8/1	
MD 12 VARIO select		
Max. pumping speed	16 m ³ /h	9.4 cfm
Ultimate vacuum, abs.	1,5 mbar	1.1 Torr
Number of cylinders/stages	8/3	
MV 10 VARIO select		
Max. pumping speed	14 m ³ /h	8.4 cfm
Ultimate vacuum, abs.	0,3 mbar	0.2 Torr
Number of cylinders/stages	8/4	
Max. inlet pressure, abs.	1,1 bar	825 Torr
Max. outlet pressure, abs.	1,1 bar	825 Torr
Max. differential pressure, abs.	1,1 bar	825 Torr
Sensor	integrated	integrated
Measuring principle	Ceramic diaphragm (aluminum oxide), capacitive, gas type independent, absolute pressure	
Accuracy of measurement	<±1 mbar/hPa/Torr, ±1 digit (after adjustment, constant temperature)	
Upper measurement limit	1080 mbar	810 Torr
Lower measurement limit	0,1 mbar	0.1 Torr
Temperature coefficient	< 0,15 mbar/hPa/K	< 0.11 Torr/K

Weights* and dimensions (l x b x h)			(US)
ME 16 VARIO select	552 mm x 260 mm x 450 mm		21.7 in x 10.2 in x 17.7 in
Weight*	28,3 kg		62.4 lb
MD 12 VARIO select	552 mm x 260 mm x 450 mm		21.7 in x 10.2 in x 17.7 in
Weight*	28,3 kg		62.4 lb
MV 10 VARIO select	552 mm x 260 mm x 450 mm		21.7 in x 10.2 in x 17.7 in
Weight*	28,3 kg		62.4 lb

* without cable

Other information	
Sensor type	VACUU-SELECT Sensor
Controller	VACUU-SELECT
Sound pressure level at 1500 rpm/62% (VARIO)	50 dBA

8.1.2 Wetted materials

Wetted materials

Component	Wetted materials
Pump	
Housing cover	Aluminium alloy
Head cover	Aluminium alloy (AlSi12)
Diaphragm clamping disc	Aluminium alloy (AlSi12)
Diaphragm	FPM
Valves	FPM
O-rings	FPM
Small flange	Stainless steel
Hose fittings	ETFE/ECTFE
Hoses	PTFE
Inlet	Aluminium alloy
Suction/pressure distributor	Aluminium alloy
Outlet	PP
Silencer	Anodized aluminum / PTFE / PTFE carbon reinforced / spring steel
VACUU-SELECT Sensor	
Vacuum sensor	Aluminum oxide ceramic, gold-coated (if applicable)
Measurement chamber	PPS
Small flange	PP
Sealing ring at the sensor	chemically resistant fluoroelastomer
Hose nozzle	PP

8.1.3 Rating plate


Data on rating plate



- ⇒ In the event of an error, make a note of the type and serial number on the rating plate.
- ⇒ When contacting our Service Department, please provide the type and serial number from the rating plate. This will allow us to provide you with specific support and advice for your device.

Diaphragm pump rating plate, general

→ Example
Cut-out showing
rating plate

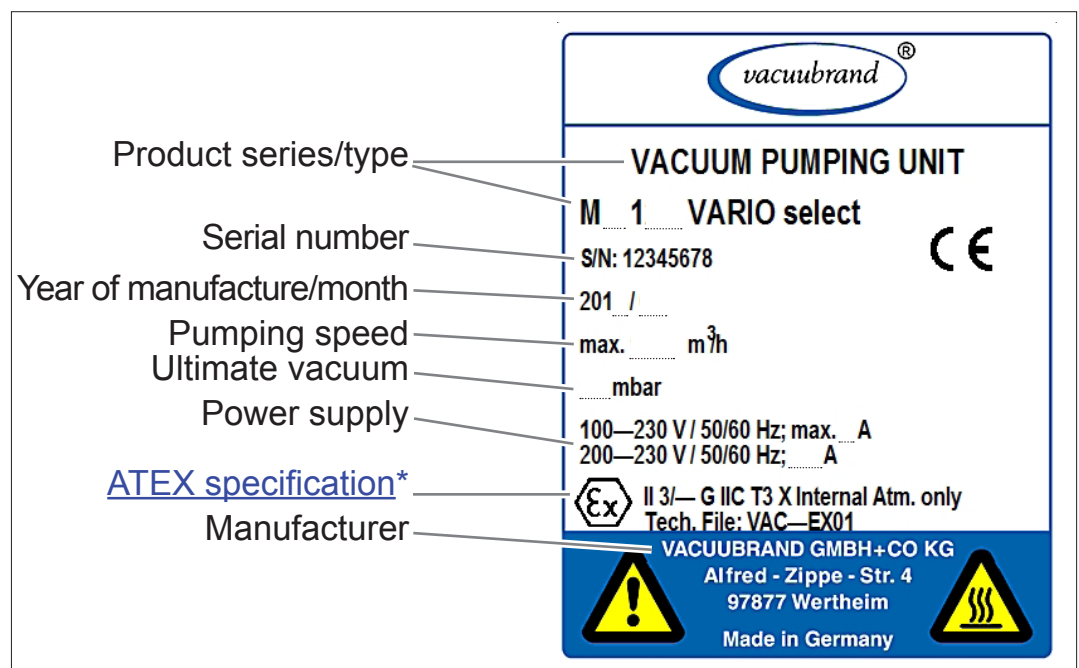

* Documentation, group and category, marking G (gas), type protection, explosion group, temperature class (additionally see : [Approval for ATEX equipment](#)).

8.2 Ordering information

Ordering information
pump series

Diaphragm pump series	Order no.*
ME 16 VARIO select	207411xx
MD 12 VARIO select	207431xx
MV 10 VARIO select	207441xx

* Order no. depends on power cable CEE, CH, UK, US, CN, IN

Ordering information
accessories

Accessories	Order no.
Separator flask AK	20699979
PTFE hose KF DN 25 (l = 1000 mm)	20686033
Hose (rubber) d _i 15 mm (length to order)	20686003
Stainless steel hose KF DN 25 (l = 1000 mm)	20673337
Coolant valve VKW-B	20674220
Venting valve VBM-B	20674217
Vacuum valve VS 25, KF DN 25	20665005
VACUU·SELECT® Sensor without venting valve	20700021
Vacuum sensor VSK 3000	20640530
Adapter small flange KF DN 25 to hose nozzle DN 15	20662808
Threaded flange KF DN 16 / 1/2"	20672101
Hose nozzle DN 15 mm / 1/2"	20642472
Adapter KF DN 25 to 2x PTFE tube DN 10/8	20667052
Adapter small flange KF DN 16 to hose nozzle 1/2"	20636004
VACUU·BUS Y adapter	20636656
Extension cable VACUU·BUS, 0.5 m	20612875
Extension cable VACUU·BUS, 2 m	20612552
Extension cable VACUU·BUS, 10 m	22618493
VACUU·BUS wall duct	20636153
DAkkS calibration with first delivery	20900214
DAkkS recalibration	20900215

Ordering information
spare parts

Spare parts		Order no.
Anti-rotation protection D17x17.5		20635113
Service kit MD 12 / MV 10 VARIO select		20696827
Service kit ME 16 VARIO select		20696819
Power cable	CEE	20612058
	CH	20676021
	CN	20635997
	IND	20635365
	UK	20612065
*Silencer G 1/2"		20642473

* *Caution: Gases containing dust, deposits and condensed solvent vapors can affect the flow of gas through the silencer. These factors or a high gas flow rate can cause excess pressure to build up, which can damage the pump bearings, diaphragms, and valves. Do not use the silencer in such circumstances.*



⇒ A full list of spare parts available can be found under
→ [VACUUBRAND > Support > Repair instructions > Diaphragm pumps.](#)

Sources of supply

Purchase original accessories and original spare parts from a subsidiary of **VACUUBRAND GMBH + CO KG** or your local distributor.

International
sales offices and
specialized trade



- ⇒ Information about our complete product range is available in the current [product catalog](#).
- ⇒ Your local distributor or **VACUUBRAND GMBH + CO KG** [sales office](#) is available to assist you with orders, questions on vacuum control and optimal accessories.

8.3 Service

Service offer and service range

Take advantage of the comprehensive range of services available from **VACUUBRAND GMBH + CO KG**.



Services in detail

- Product consultation and practical solutions
- Fast delivery of spare parts and accessories
- Professional maintenance
- Immediate repairs processing
- On-site service (on request)
- [Calibration](#) (DAkkS-accredited)
- With Health and Safety Clearance form: Return, disposal.


⇒ Visit our website for further information: www.vacuubrand.com.

Service handling

Meet the terms of service

1. Contact your local distributor or our Service Department.
2. Request an RMA no. for your order.
3. Clean the product thoroughly or if necessary, decontaminate it professionally.
4. Fill out the [Health and Safety Clearance form](#) in full.
5. Return your product, including:
 - RMA no. and description of the error
 - Repair or service order,
 - Health and Safety Clearance form
 - Attach everything to the outside of the package

Return (reshipment)

⇒  Reduce downtime, speed up processing. Please keep the required data and documents ready when contacting our Service Department.

- ▶ Your order can be quickly and easily processed.
- ▶ Hazards can be prevented.
- ▶ A brief description and/or photos will help locate the source of the error.

8.4 Index

Index	<p>A</p> <p>Abbreviations 9</p> <p>Accessories 75</p> <p>Action steps 9</p> <p>Additional symbols 8</p> <p>Assemble silencer 30</p> <p>ATEX equipment category 21</p> <p>ATEX equipment labeling 21</p> <p>B</p> <p>Blockage protection 20</p> <p>C</p> <p>Clean the surfaces 50</p> <p>Copyright © 5</p> <p>D</p> <p>Device and housing parts, assemble 68</p> <p>Device and housing parts, disassemble 52</p> <p>Diaphragm pump stages 24</p> <p>Disassemble the front housing section 52</p> <p>Display conventions 7</p> <p>Display of operating steps 9</p> <p>E</p> <p>EC Declaration of Conformity 80</p> <p>Electrical connection 35</p> <p>Eliminate sources of danger 17</p> <p>Error – Cause – Remedy ... 43, 44, 45</p> <p>Exhaust gas connection 33</p> <p>Explanation of safety symbols 8</p> <p>Exploded-drawing pump head 55</p> <p>Exploded-drawing suction/pressure distributor 66</p> <p>F</p> <p>Foreseeable misuse 12, 13</p> <p>G</p> <p>Goods arrival 25</p> <p>H</p> <p>Handling instructions 9</p> <p>Hazards during venting 18</p> <p>Hot surfaces 19</p> <p>I</p> <p>Icons 8</p> <p>Improper use 12</p> <p>Incorrect measurements 17</p> <p>Insert valves and O-rings 61</p> <p>Instruction modules 6</p> <p>Intended use 11</p> <p>L</p> <p>Limitation of use 27</p>	<p>M</p> <p>Maintain minimum distance 19</p> <p>Maintain the top, right and left pump head pairs 62</p> <p>Maintenance 47</p> <p>Maintenance intervals 48</p> <p>Maintenance of suction/pressure distributor 65</p> <p>Mandatory sign 8</p> <p>Manual structure 6</p> <p>Measurement chamber 73</p> <p>Minimum distances 27</p> <p>Mount the base (option) 28</p> <p>Mount the side panel 63</p> <p>O</p> <p>Operating elements Vacuum controller 38</p> <p>Operating hours until maintenance 48</p> <p>Operating panel 10</p> <p>Operator obligations 13</p> <p>Ordering information 75</p> <p>Ordering information pump series .. 75</p> <p>Overheating 19</p> <p>Overheating protection 20</p> <p>P</p> <p>Peripheral devices and ATEX 21</p> <p>Personnel obligations 13</p> <p>Personnel qualification 14</p> <p>Power cable 35</p> <p>Pressure relief valve 65</p> <p>Pressure relief valve + O-ring, maintenance 65</p> <p>Pressure relief valve + O-ring, replacement 66</p> <p>Prevent blockages in the outlet line 17</p> <p>Prevent condensate return 17</p> <p>Prevent ignition sources 21</p> <p>Process screen (main screen) 38</p> <p>Product-specific terms 10</p> <p>Prohibition sign 8</p> <p>Proper disposal 22</p> <p>Protective clothing 15</p> <p>Pump electrical connection 35</p> <p>Q</p> <p>Qualification description 14</p> <p>Quality standard 15</p> <p>R</p> <p>Rating plate 74</p> <p>Recommended aids for cleaning and maintenance 48</p> <p>Remove the left side panel 53</p> <p>Remove the right side panel 54</p> <p>Reset thermal fault 44</p> <p>Residual energies 18</p> <p>Responsibility Assignment Matrix .. 14</p>
-------	--	---

Index

- Responsibility matrix and areas of competence 14
- Return (reshipment) 77
- Return to use after maintenance . . . 69

S

- Safety information 11
- Safety information for vacuum equipment 11
- Safety precautions 15
- Sales offices 76
- Service handling 77
- Service range 77
- Silencer 29
- Silencer at the outlet EX 29
- Sources of supply 76
- Spare parts 76
- Start 39
- Stop 39
- Storage, prepared 69
- Switch pump off 40
- Switch pump on 37
- Symbols 8

T

- Target groups 14
- Technical data 70, 71
- Technical information 70
- Technical support 43

U

- Unpacking 25
- User information 5

V

- VACUU·BUS® 10
- VACUU·BUS® connector 10
- Vacuum connection 31
- Vacuum controller operation 37
- Valves, pump head 55
- Venting 39

W

- Warning sign 8
- Warning signs and labels 20
- Wetted materials 73

8.5 EC Declaration of Conformity

**EU-Konformitätserklärung
EC Declaration of Conformity
Déclaration CE de conformité**



Hersteller / Manufacturer / Fabricant:

VACUUBRAND GMBH + CO KG · Alfred-Zippe-Str. 4 · 97877 Wertheim · Germany

Hiermit erklärt der Hersteller, dass das Gerät konform ist mit den Bestimmungen der Richtlinien:

Hereby the manufacturer declares that the device is in conformity with the directives:

Par la présente, le fabricant déclare, que le dispositif est conforme aux directives:

2006/42/EG (M-RL), 2014/30/EU (EMV-RL), 2014/34/EU (ATEX-RL),
2011/65/EU (RoHS-2)

Membranpumpen / Diaphragm pumps / Pompes à membrane

Typ / Type / Type: **ME 16 VARIO select, MD 12 VARIO select, MV 10 VARIO select**

Artikelnummer / Order number / Numéro d'article: **20741150, 20743150,
20744150**

Seriennummer / Serial number / Numéro de série: Siehe Typenschild / See rating plate / Voir plaque signalétique

Angewandte harmonisierte Normen / Harmonized standards applied / Normes harmonisées utilisées: DIN EN ISO 12100:2011, DIN EN 1012-2:2011, IEC 61010-1:2010 (Ed. 3), DIN EN 61010-1:2011, DIN EN 61326-1:2013, DIN EN 1127-1:2011, DIN EN ISO 80079-36:2016, DIN EN 50581:2013

Bevollmächtigter für die Zusammenstellung der technischen Unterlagen / Person authorised to compile the technical file / Personne autorisée à constituer le dossier technique: Dr. J. Dirscherl · VACUUBRAND GMBH + CO KG · Germany

Ort, Datum / place, date / lieu, date: Wertheim, 02.08.2018



(Dr. F. Gitmans)

*Geschäftsführer / Managing Director /
Gérant*

ppa.



(Dr. J. Dirscherl)

*Technischer Leiter / Technical
Director / Directeur technique*

VACUUBRAND GMBH + CO KG

Alfred-Zippe-Str. 4
97877 Wertheim

Tel.: +49 9342 808-0

Fax: +49 9342 808-5555

E-Mail: info@vacuubrand.com

Web: www.vacuubrand.com



Technology for Vacuum Systems

Manufacturer:

VACUUBRAND GMBH + CO KG
Alfred-Zippe-Str. 4
97877 Wertheim
GERMANY

Phone:

- Head office +49 9342 808-0
- Sales +49 9342 808-5550
- Service +49 9342 808-5660

Fax: +49 9342 808-5555

Email: info@vacuubrand.com

Web: www.vacuubrand.com